

Legacy™ System Catalog



Legacy2



Legacy1



Legacy2



Legacy3






Table of Contents

	Page #	
	1-5	Implant Design, Surface and Dimensions
	6-11	Legacy1, Legacy2 and Legacy3 Implants 
	12-13	Labeling and Fixture-Mount Packaging
	14	Procedures for Contoured Abutments 
	15	Procedures for Titanium FMT / Temporary Abutment
	16-17	Selection Chart for Legacy Systems 
	18	Cover Screws, Healing Collars, Transfers & Analogs
	19	Straight Abutment & Ancillary Components 
	20	Straight & Angled Contoured, Titanium & Zirconia Abutments
	21	Gold, Engaging and Non-Engaging Abutments 
	22	Selection Chart for Overdenture & Screw-Retained Restorations
	23	Screw-Receiving, Ball and GPS Abutments 
	24	Overdenture and Screw-Receiving Abutment Ancillary Products
	25-27	Instruments, Surgical Drills and Surgical Sequence 
	28-29	Surgical Tray



Mini-threads with Double-lead Threads for Faster Insertion

The bodies of the Legacy1, Legacy2 and Legacy3 implants are tapered to facilitate insertion in an undersized socket (stop at intermediate drill) for expansion in soft bone. This increases initial stability (G. Niznick: *Clinical Oral Health 2000 "Achieving Success in Soft Bone"*). The SBM textured surface extends over the entire endosseous portion of the Legacy system implants. Mini-threads start 1mm below the crest of the bone and extend for 2.0-2.5mm before transitioning into double-lead threads to the apex of the implant. The mini-threads are half the depth and half the distance between adjacent peaks compared to the double-lead threads that extend to the apex of the implant. Mini-threads reduce stress in the critical crestal bone region while double-lead threads reduce the number of turns required to seat the implant. Vertical cutting grooves extend up the implant from the apex for self-tapping insertion.

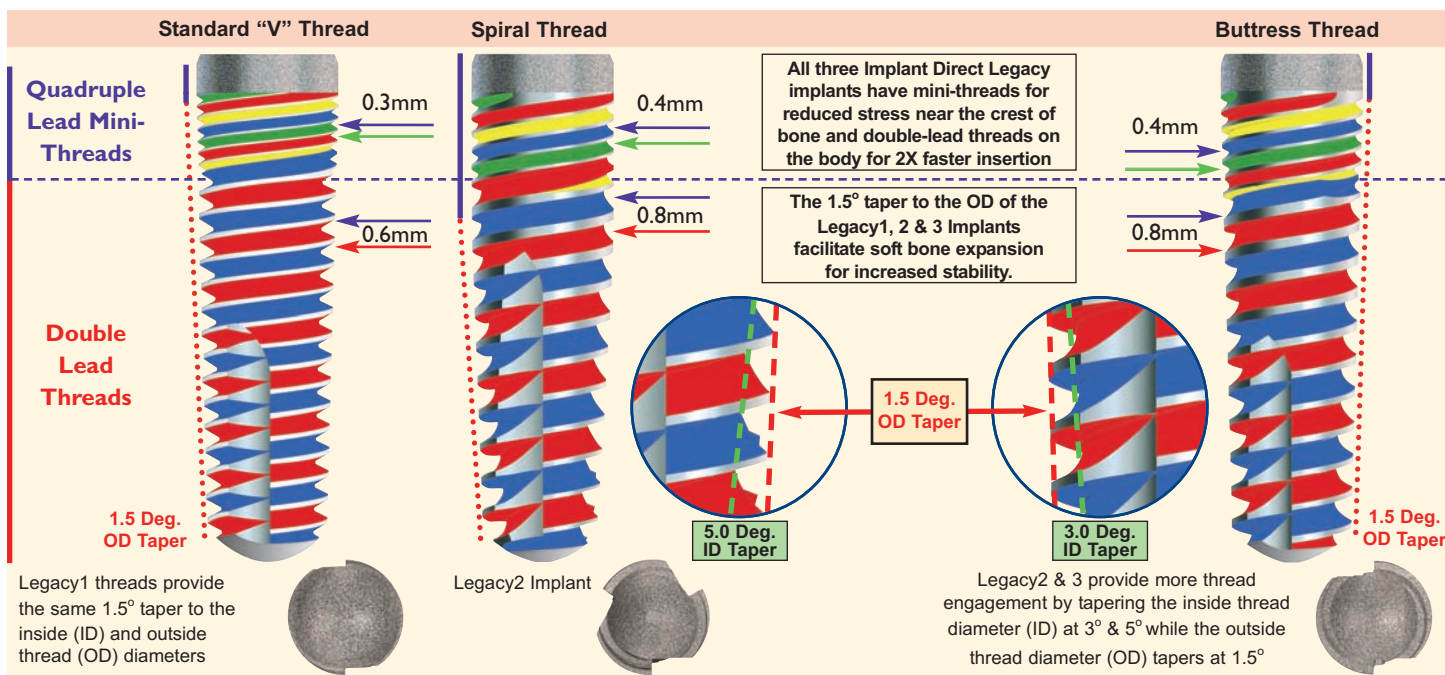


SBM and HA Coated Surface Options

SBM surface is created by using a **S**oluble **B**last **M**edia of Hydroxylapatite (HA) crystals to create a medium-rough texture. Implant Direct uses the same vendor for its SBM and HA coated surface treatments as used by Paragon Implant Company from 1992-2000 when these studies were conducted.

SEM's of Implant Direct's Surfaces


Mini-threads textured by blasting with **SBM™**: Blasted with Soluble HA Particles
HA Coating: High Crystalline 40u layer

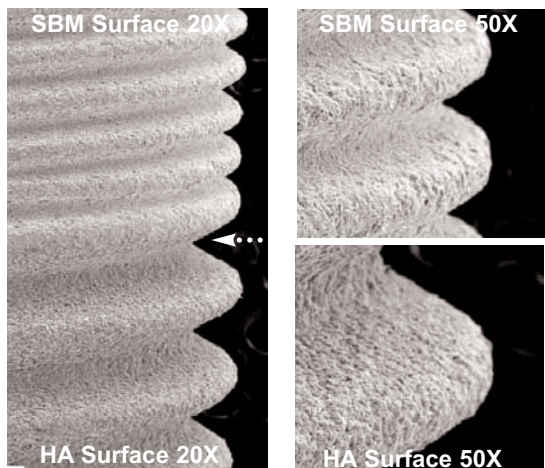


Animal Studies Bone Response to SBM Surface



Brazilian Study* Shows SBM Provides Superior Bone Contact
 Mean Percentage of Bone Contact


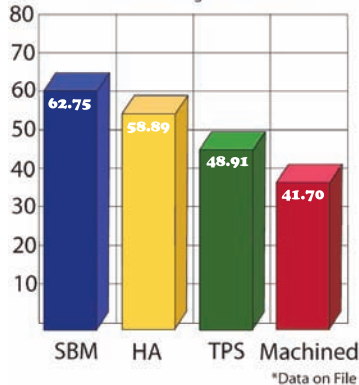


LSU Study* Shows SBM Provides Superior Bone Contact


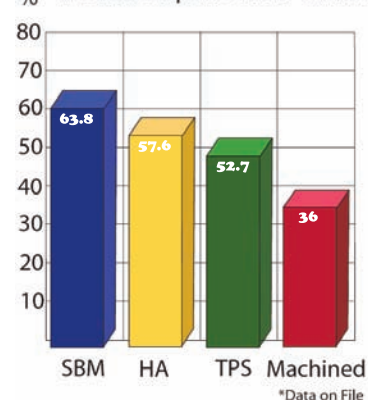






Legacy I Implant Dimensions

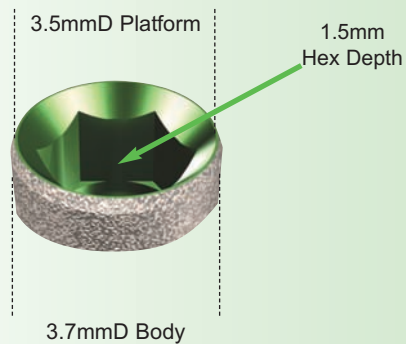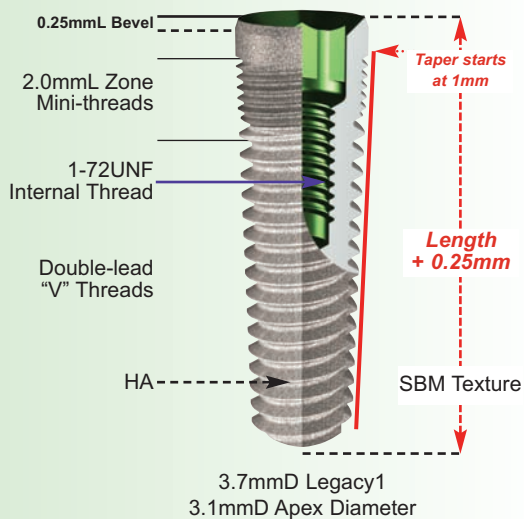

3.5mmD Platform with a 2.5mmD Internal Hex

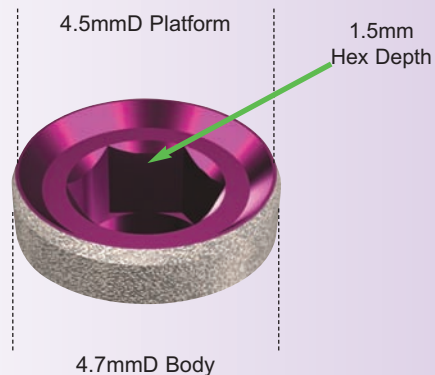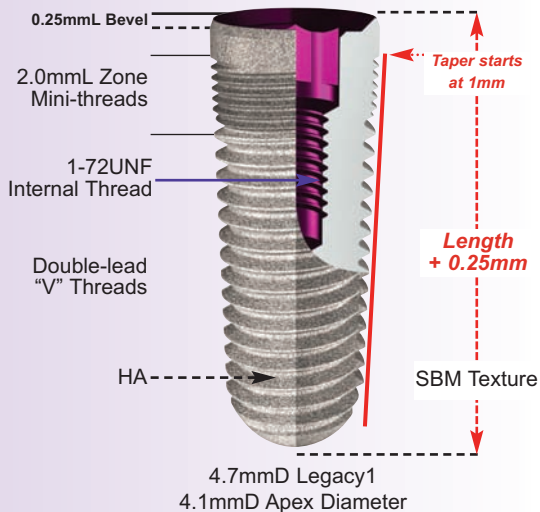









4.5mmD Platform with a 2.5mmD Internal Hex

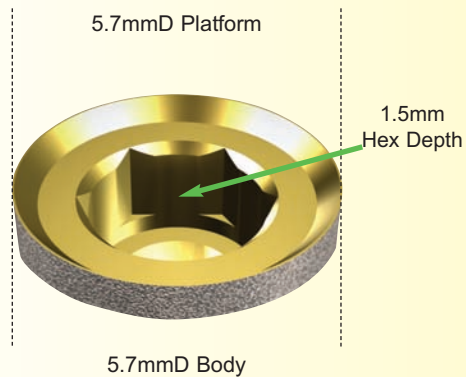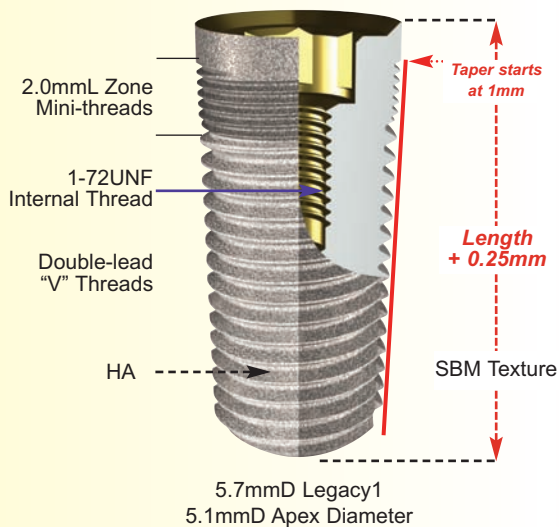

5.7mmD Platform with a 3.0mmD Internal Hex





Legacy2 Implant Dimensions

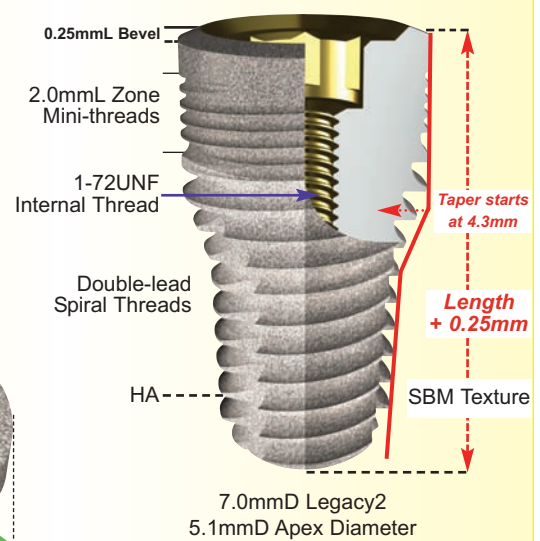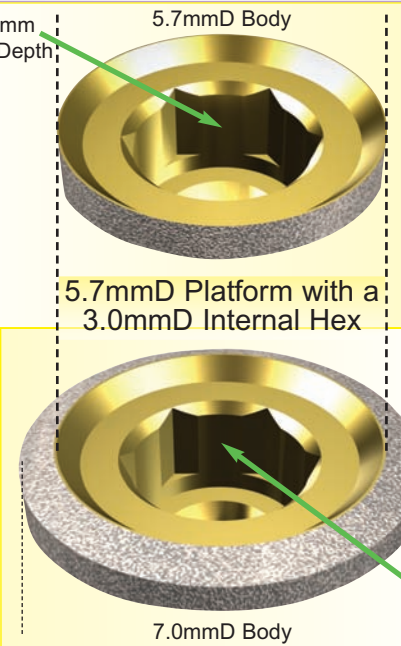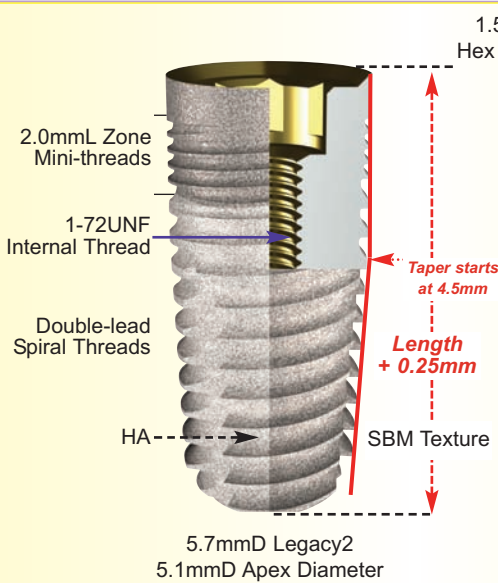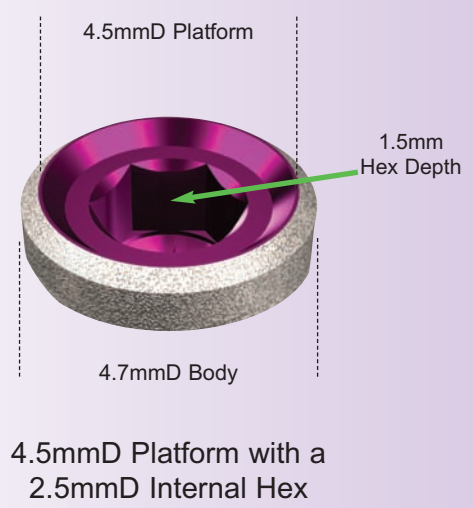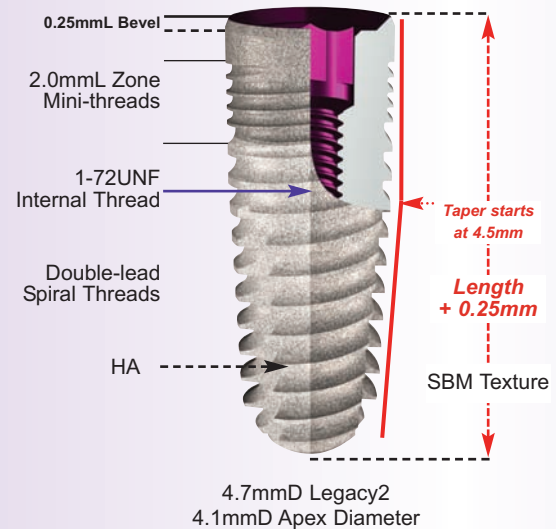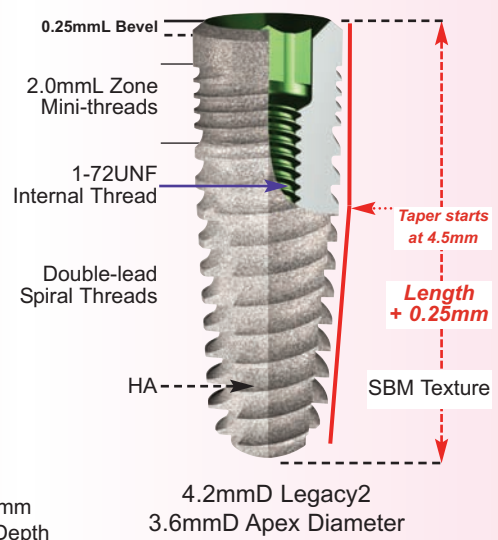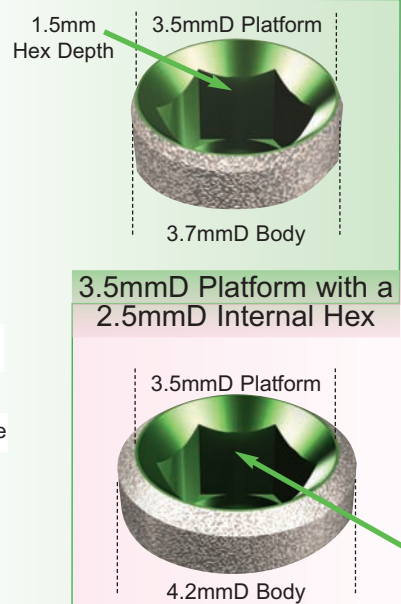
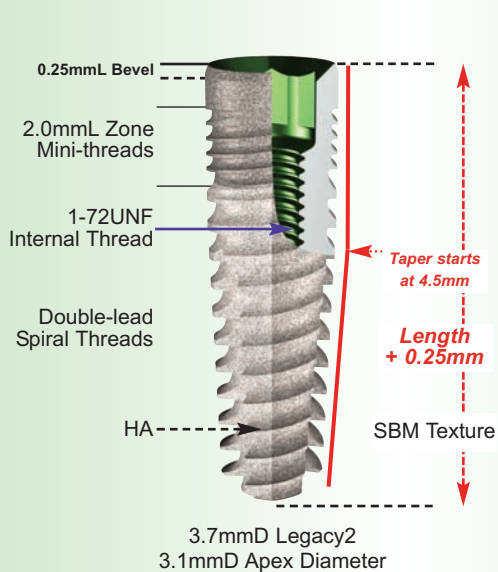






Legacy3 Implant Dimensions

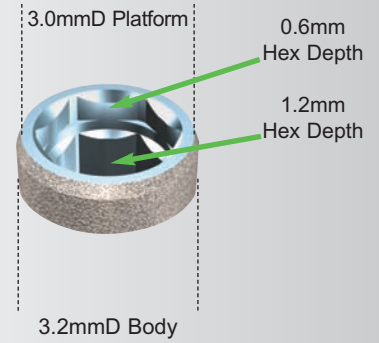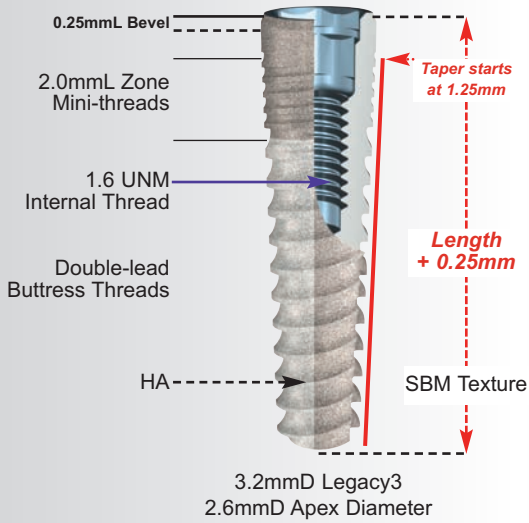

3.0mmD Platform with a Double Internal Hex
2.25mmD (upper) and 2.0mmD (lower)

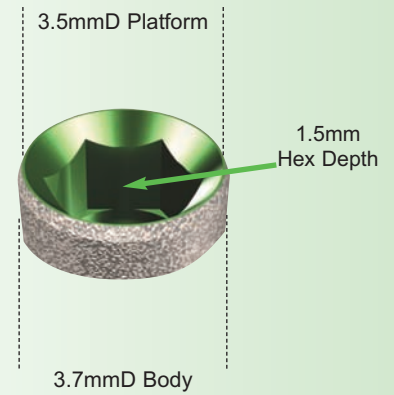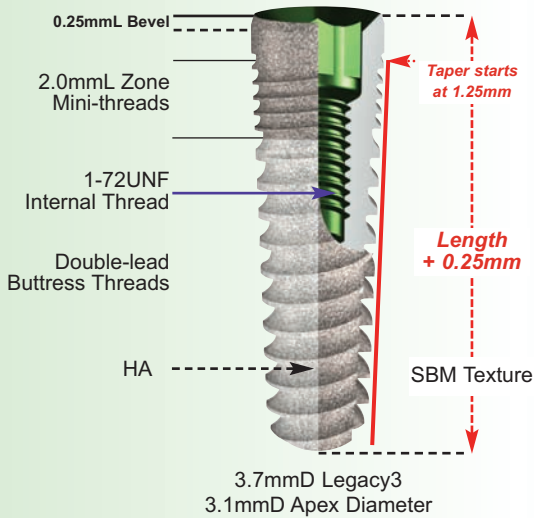



3.5mmD Platform with a 2.5mmD Internal Hex

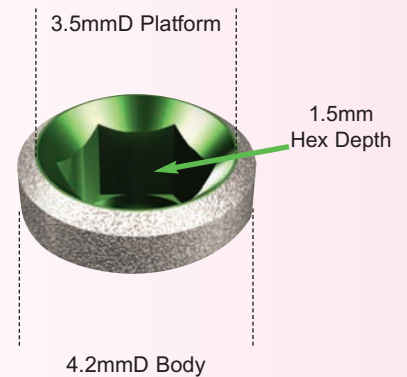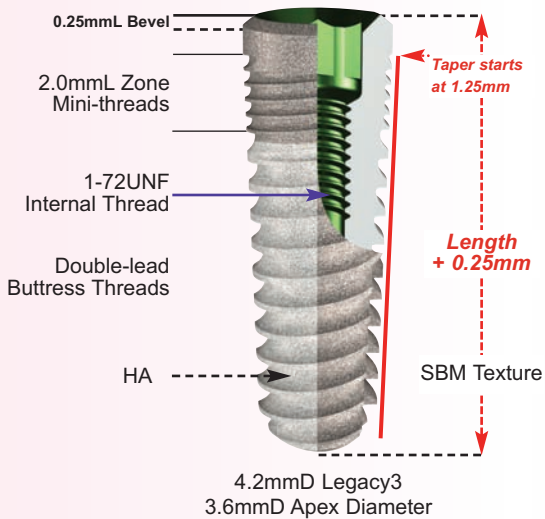



3.5mmD Platform with a 2.5mmD Internal Hex





Legacy3 Implant Dimensions

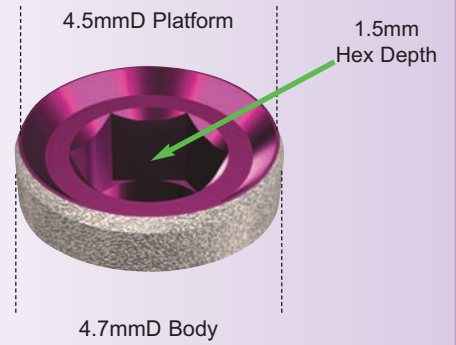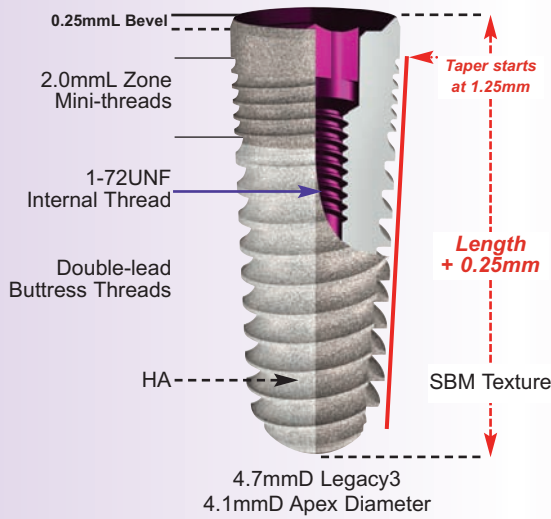

4.5mmD Platform with a 2.5mmD Internal Hex

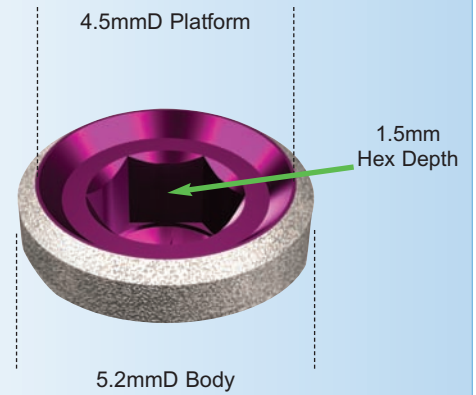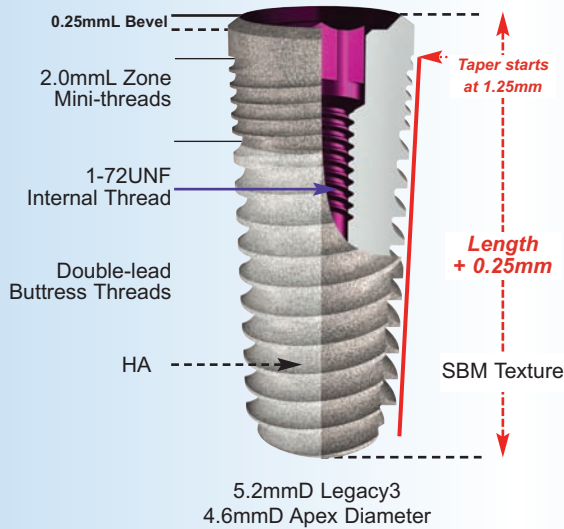

4.5mmD Platform with a 2.5mmD Internal Hex

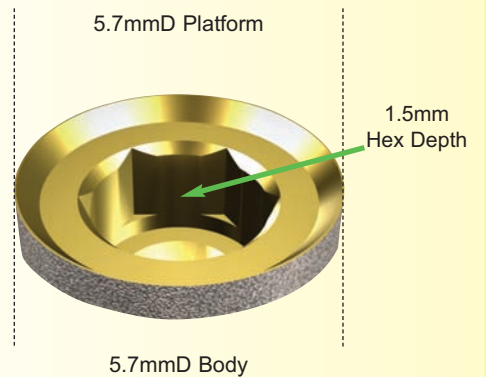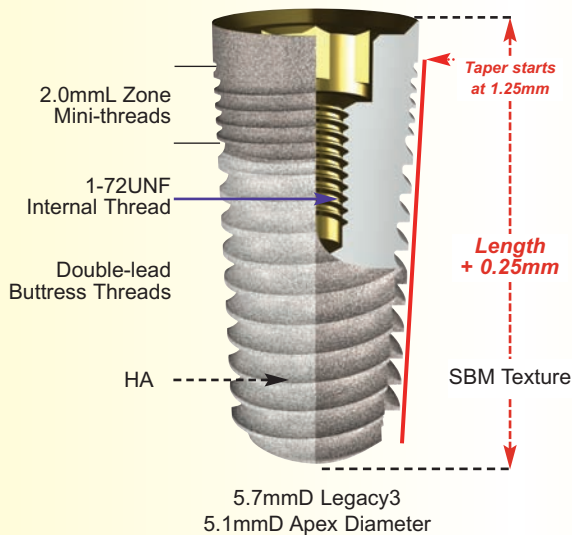

5.7mmD Platform with a 3.0mmD Internal Hex


Legacy I Implant System Features and Benefits

- Internal connection with 1.5mm long hex and internal bevel
- Tapered body allows insertion using Tapered Screw-Vent® drills (except the 5.7mmD implant uses 5.4/4.8mmD Step Drill)
- Surface options: SBM (Soluble Blast Media) textured or HA-coated
- Double-lead “V” threads for 2X faster insertion and quadruple-lead threads near top for increased thread engagement and reduced stress
- 2mm extender retained by cover screw for use as a healing collar for one stage surgery
- Two self-tapping grooves, extending halfway up implant






Legacy I Packaging

4.5mmD Platform
804713
13mm Length

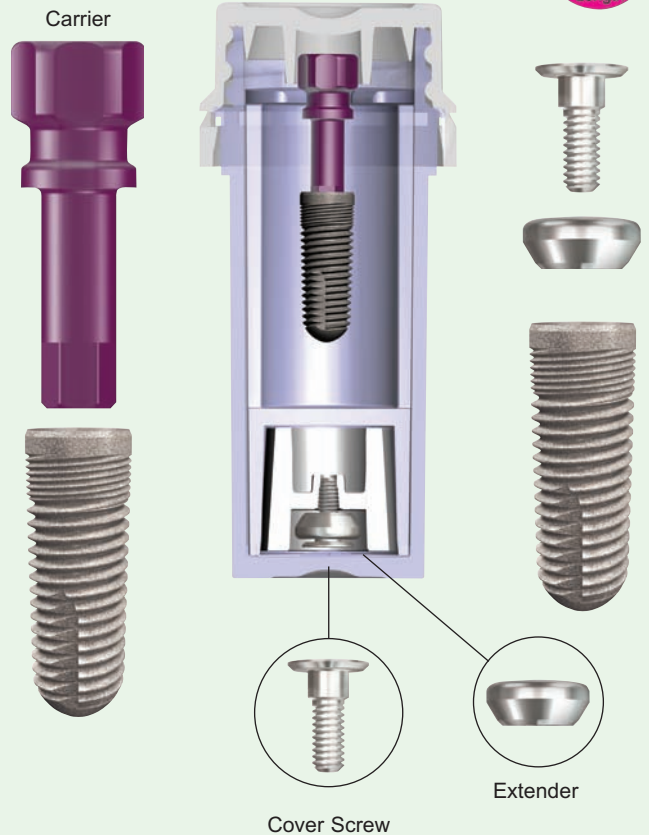

Placing Implant with Plastic Carrier

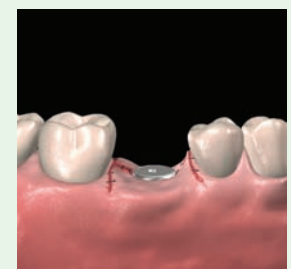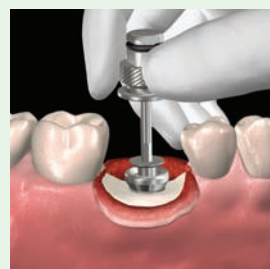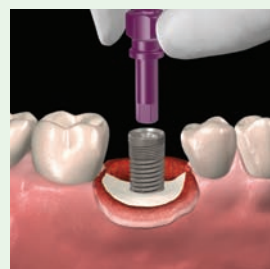



Legacy I: Two-Piece Implant for Two-Stage Surgery

Legacy I

SBM and HA Surface Implants Includes Cover Screw & Extender

3.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
803708	\$125.00	813708	\$150.00	3.7mmD	8mmL	3.5mmD
803710	\$125.00	813710	\$150.00	3.7mmD	10mmL	3.5mmD
803711	\$125.00	813711	\$150.00	3.7mmD	11.5mmL	3.5mmD
803713	\$125.00	813713	\$150.00	3.7mmD	13mmL	3.5mmD
803716	\$125.00	813716	\$150.00	3.7mmD	16mmL	3.5mmD



4.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
804708	\$125.00	814708	\$150.00	4.7mmD	8mmL	4.5mmD
804710	\$125.00	814710	\$150.00	4.7mmD	10mmL	4.5mmD
804711	\$125.00	814711	\$150.00	4.7mmD	11.5mmL	4.5mmD
804713	\$125.00	814713	\$150.00	4.7mmD	13mmL	4.5mmD
804716	\$125.00	814716	\$150.00	4.7mmD	16mmL	4.5mmD



5.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
805708	\$125.00	815708	\$150.00	5.7mmD	8mmL	5.7mmD
805710	\$125.00	815710	\$150.00	5.7mmD	10mmL	5.7mmD
805711	\$125.00	815711	\$150.00	5.7mmD	11.5mmL	5.7mmD
805713	\$125.00	815713	\$150.00	5.7mmD	13mmL	5.7mmD
805716	\$125.00	815716	\$150.00	5.7mmD	16mmL	5.7mmD




Legacy2 Implant System Features and Benefits


- Internal connection with 1.5mm long hex and internal lead-in bevel
- Body is straight for top 1/3rd with bottom 2/3rds tapered for bone expansion. Surgically compatible with Tapered Screw-Vent® drills (except the 5.7mmD Implant uses 5.4/4.8mmD Step Drill)
- Surface options: SBM (Soluble Blast Media) textured or HA-coated
- Double-lead spiral threads for improved insertion and quadruple-lead threads near top for increased thread engagement and reduced stress
- 2mm extender retained by cover screw for use as a healing collar for one stage surgery
- Color-coded carrier usable as a closed-tray transfer and temporary abutment
- Three aggressive self-tapping grooves extending 2/3rds length up the implant

Legacy2 Packaging


3.5mmD
Platform
823713
13mm
Length

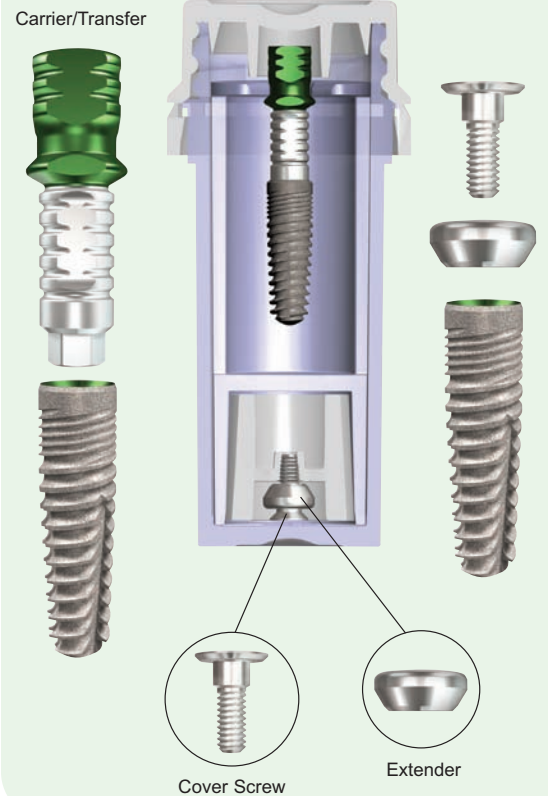


3.5mmD
Platform

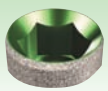


3.5mmD
Platform

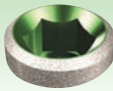

4.5mmD
Platform

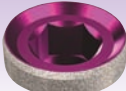

5.7mmD
Platform

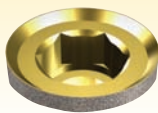

5.7mmD
Platform



3.7mmD
Implant




4.2mmD
Implant



4.7mmD
Implant



5.7mmD
Implant

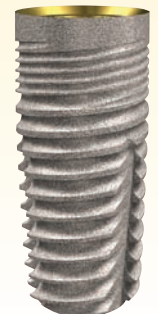

7.0mmD
Implant




Placing Carrier/Transfers & Making Impressions




Color-coded titanium fixture-mount carrier



Align flat to labial for indexing of hex

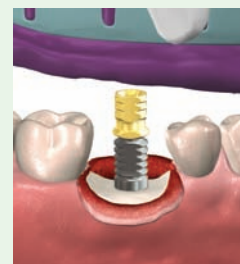

Use fixture-mount as closed tray transfer

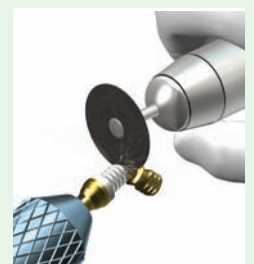

Shorten fixture-mount for temporary abutment



Legacy2: Two-Piece Implant for Two-Stage Surgery

Legacy2

SBM and HA Surface Implants Includes Cover Screw, Extender & Carrier/Transfer

3.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
823708	\$150.00	833708	\$175.00	3.7mmD	8mmL	3.5mmD
823710	\$150.00	833710	\$175.00	3.7mmD	10mmL	3.5mmD
823711	\$150.00	833711	\$175.00	3.7mmD	11.5mmL	3.5mmD
823713	\$150.00	833713	\$175.00	3.7mmD	13mmL	3.5mmD
823716	\$150.00	833716	\$175.00	3.7mmD	16mmL	3.5mmD



4.2mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
824208	\$150.00	834208	\$175.00	4.2mmD	8mmL	3.5mmD
824210	\$150.00	834210	\$175.00	4.2mmD	10mmL	3.5mmD
824211	\$150.00	834211	\$175.00	4.2mmD	11.5mmL	3.5mmD
824213	\$150.00	834213	\$175.00	4.2mmD	13mmL	3.5mmD
824216	\$150.00	834216	\$175.00	4.2mmD	16mmL	3.5mmD



4.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
824708	\$150.00	834708	\$175.00	4.7mmD	8mmL	4.5mmD
824710	\$150.00	834710	\$175.00	4.7mmD	10mmL	4.5mmD
824711	\$150.00	834711	\$175.00	4.7mmD	11.5mmL	4.5mmD
824713	\$150.00	834713	\$175.00	4.7mmD	13mmL	4.5mmD
824716	\$150.00	834716	\$175.00	4.7mmD	16mmL	4.5mmD



5.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
825708	\$150.00	835708	\$175.00	5.7mmD	8mmL	5.7mmD
825710	\$150.00	835710	\$175.00	5.7mmD	10mmL	5.7mmD
825711	\$150.00	835711	\$175.00	5.7mmD	11.5mmL	5.7mmD
825713	\$150.00	835713	\$175.00	5.7mmD	13mmL	5.7mmD
825716	\$150.00	835716	\$175.00	5.7mmD	16mmL	5.7mmD



7.0mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
827008	\$150.00	837008	\$175.00	7.0mmD	8mmL	5.7mmD
827010	\$150.00	837010	\$175.00	7.0mmD	10mmL	5.7mmD
827011	\$150.00	837011	\$175.00	7.0mmD	11.5mmL	5.7mmD
827013	\$150.00	837013	\$175.00	7.0mmD	13mmL	5.7mmD



Legacy3 Implant System Features and Benefits

- Internal connection with 1.5mm long hex and internal lead-in bevel
- Tapered body allows insertion using Tapered Screw-Vent® drills for 3.7mmD, 4.2mmD and 4.7mmD implants only
- Surface options: SBM (Soluble Blast Media) textured or HA-coated
- Double-lead buttress threads for improved insertion and initial stability, quadruple-lead threads near top for increased thread engagement and reduced stress
- 2mm extender retained by cover screw for use as a healing collar or for one or two-stage surgery
- Color-coded carrier usable as a closed-tray transfer and/or final full-contour abutment
- Two self-tapping grooves, extending halfway up implant

Legacy3 Packaging

4.5mmD
Platform
854713
13mm
Length

Carrier/Transfer



Cover Screw




Extender


3.0mmD
Platform



3.5mmD
Platform

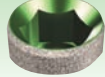

3.5mmD
Platform

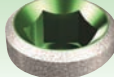

4.5mmD
Platform



4.5mmD
Platform



5.7mmD
Platform



3.2mmD
Implant



3.7mmD
Implant



4.2mmD
Implant




4.7mmD
Implant



5.2mmD
Implant



5.7mmD
Implant

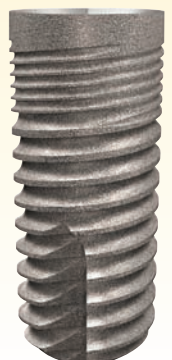



Legacy3: Two-Piece Implant for Two-Stage Surgery

Legacy3

SBM and HA Surface Implants

Includes Cover Screw, Extender & Carrier/Transfer
*Extender NOT supplied with 6mmL implant

3.2mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
853208	\$175.00	863208	\$200.00	3.2mmD	8mmL	3.0mmD
853210	\$175.00	863210	\$200.00	3.2mmD	10mmL	3.0mmD
853211	\$175.00	863211	\$200.00	3.2mmD	11.5mmL	3.0mmD
853213	\$175.00	863213	\$200.00	3.2mmD	13mmL	3.0mmD
853216	\$175.00	863216	\$200.00	3.2mmD	16mmL	3.0mmD




3.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
853706	\$175.00	863706	\$200.00	3.7mmD	6mmL	3.5mmD
853708	\$175.00	863708	\$200.00	3.7mmD	8mmL	3.5mmD
853710	\$175.00	863710	\$200.00	3.7mmD	10mmL	3.5mmD
853711	\$175.00	863711	\$200.00	3.7mmD	11.5mmL	3.5mmD
853713	\$175.00	863713	\$200.00	3.7mmD	13mmL	3.5mmD
853716	\$175.00	863716	\$200.00	3.7mmD	16mmL	3.5mmD



4.2mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
854206	\$175.00	864206	\$200.00	4.2mmD	6mmL	3.5mmD
854208	\$175.00	864208	\$200.00	4.2mmD	8mmL	3.5mmD
854210	\$175.00	864210	\$200.00	4.2mmD	10mmL	3.5mmD
854211	\$175.00	864211	\$200.00	4.2mmD	11.5mmL	3.5mmD
854213	\$175.00	864213	\$200.00	4.2mmD	13mmL	3.5mmD
854216	\$175.00	864216	\$200.00	4.2mmD	16mmL	3.5mmD

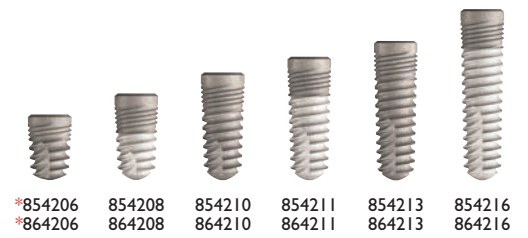

4.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
854706	\$175.00	864706	\$200.00	4.7mmD	6mmL	4.5mmD
854708	\$175.00	864708	\$200.00	4.7mmD	8mmL	4.5mmD
854710	\$175.00	864710	\$200.00	4.7mmD	10mmL	4.5mmD
854711	\$175.00	864711	\$200.00	4.7mmD	11.5mmL	4.5mmD
854713	\$175.00	864713	\$200.00	4.7mmD	13mmL	4.5mmD
854716	\$175.00	864716	\$200.00	4.7mmD	16mmL	4.5mmD




5.2mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
855206	\$175.00	865206	\$200.00	5.2mmD	6mmL	4.5mmD
855208	\$175.00	865208	\$200.00	5.2mmD	8mmL	4.5mmD
855210	\$175.00	865210	\$200.00	5.2mmD	10mmL	4.5mmD
855211	\$175.00	865211	\$200.00	5.2mmD	11.5mmL	4.5mmD
855213	\$175.00	865213	\$200.00	5.2mmD	13mmL	4.5mmD
855216	\$175.00	865216	\$200.00	5.2mmD	16mmL	4.5mmD




5.7mmD

SBM Code	Price	HA Code	Price	Diameter	Length	Platform
855706	\$175.00	865706	\$200.00	5.7mmD	6mmL	5.7mmD
855708	\$175.00	865708	\$200.00	5.7mmD	8mmL	5.7mmD
855710	\$175.00	865710	\$200.00	5.7mmD	10mmL	5.7mmD
855711	\$175.00	865711	\$200.00	5.7mmD	11.5mmL	5.7mmD
855713	\$175.00	865713	\$200.00	5.7mmD	13mmL	5.7mmD
855716	\$175.00	865716	\$200.00	5.7mmD	16mmL	5.7mmD

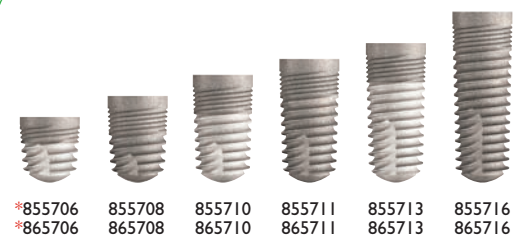




Legacy *Implant* Vial and Cap Label

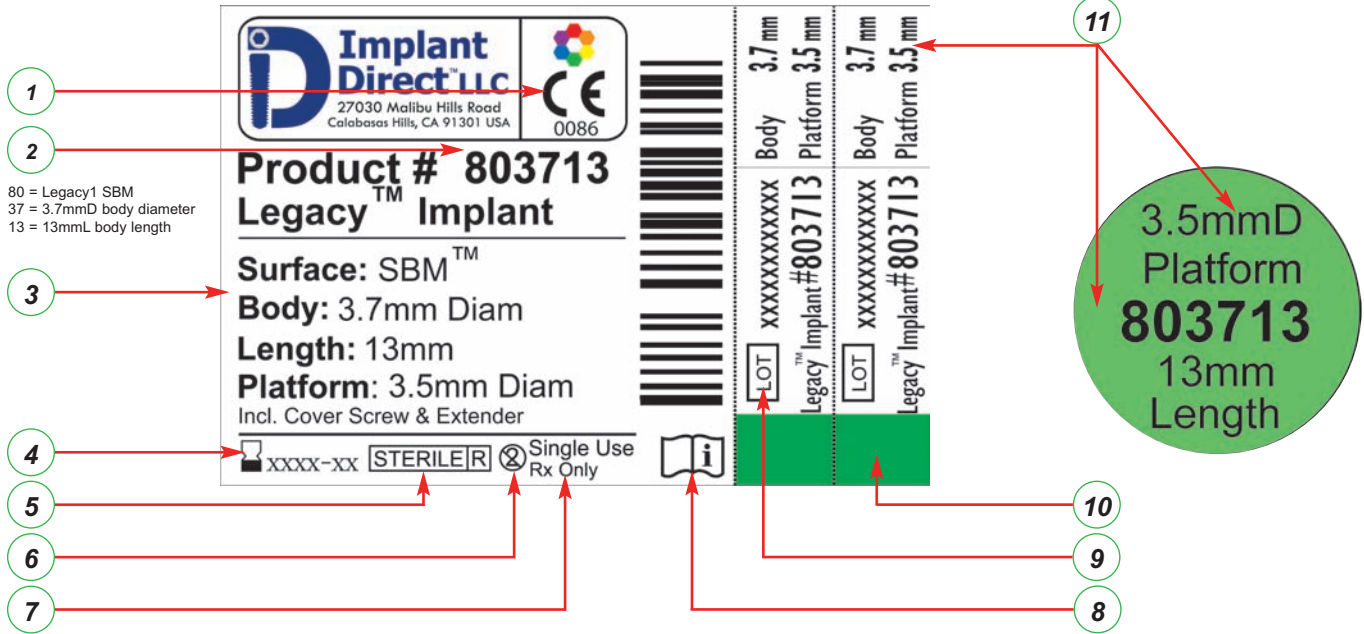







Legacy *Prosthetic* Vial and Cap Label

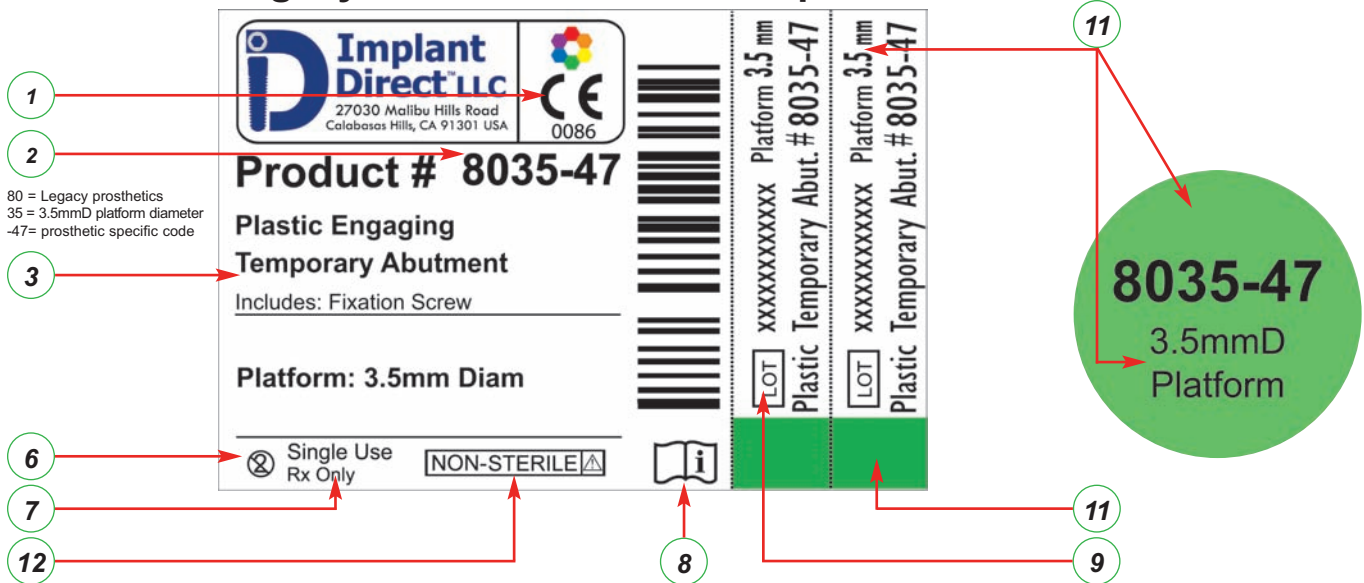

	1) Implant direct products carry the CE mark and fulfill the requirements of the Medical Device Directive 93/42/EEC
	2) Product code or reference number
	3) Product description
	4) Use before expiration date (XXXX-XX)
	5) Gamma sterilization symbol
	6) Single-use only symbol
	7) Caution: US law restricts the sale of these devices by or on the order of a dentist or physician
	8) Refer to Instructions For Use (IFU)
	9) Lot number (XXXXXXXXXX)
	10) Implant body diameter color reference on surgical tray
	11) Implant prosthetic platform diameter (written and color-coded)
	12) Non-sterile



The implant is mounted on a color-coded Carrier/Transfer.
A cover screw and extender are included in each vial.

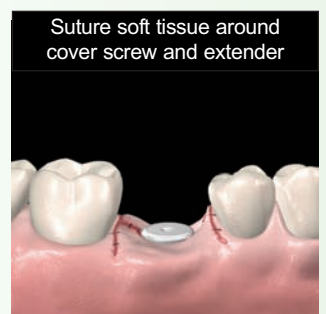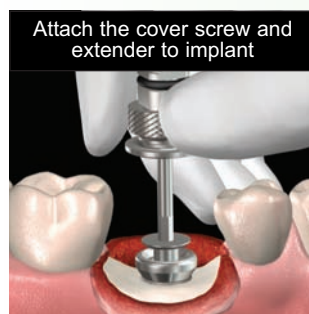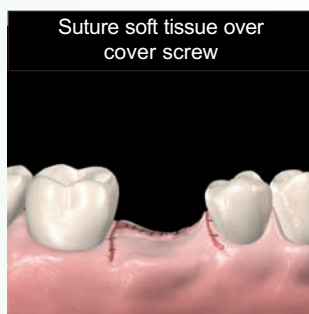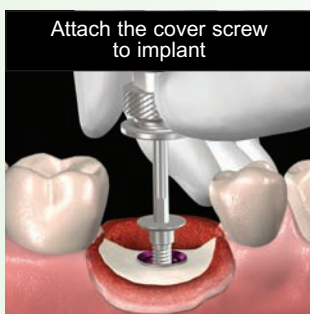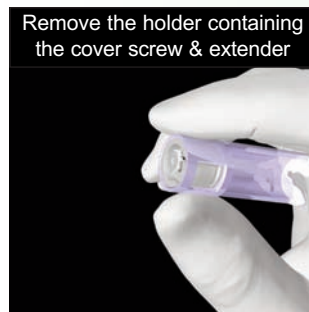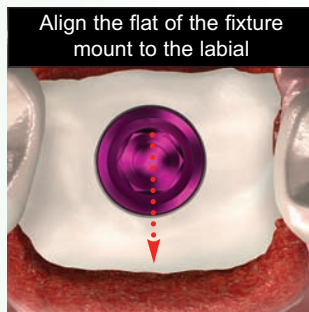



Orienting the 15° Angled Abutment



Fig. 1
Fixture Mount/
Transfer

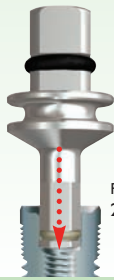


Fig. 2
2.5mmD Hex Tool


The flat side of the transfer [Fig. 1] or Hexagon Tool [Fig. 2] and the angle of the 15° abutment [Fig. 3] are manufactured to align with the flat of the implant's internal hexagon.

Fig. 3 The two-piece, 15° Angled Abutment consists of an angled body and a fixation screw




To establish the optimum angle of the abutment, index the flat surface of the Fixture Mount/Transfer either toward or away from the desired direction of abutment angle at the time of implant placement.

Option 1 : Attaching an Angled Abutment and **chairside procedures** for preparation

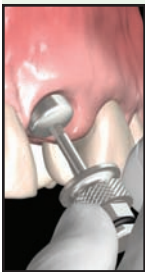


Remove the Healing Collar with a 1.25mmD Hex Tool

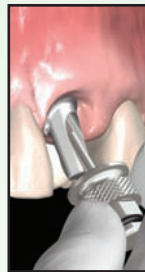

Attach the abutment with the 1.25mmD Hex Tool



Mark the required modifications to achieve appropriate clearance

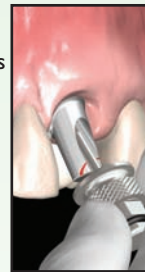


Detach the abutment with the 1.25mmD Hex Tool




Attach the abutment to the corresponding Implant Analog and reduce the abutment to achieve appropriate clearance


Option 2 : Stage-One transfer followed by **laboratory procedures** for abutment preparation and provisionalization

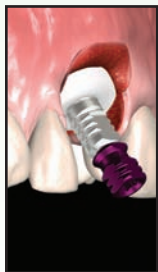

During implant placement, index the flat of the Fixture-Mount

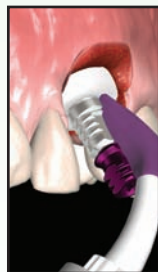


With Fixture-Mount in place, make a conventional transfer impression

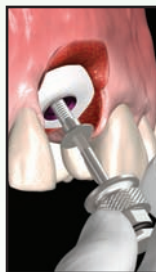


Attach the Cover Screw and submerge implant




The Fixture-Mount is retained within the impression material when it is removed

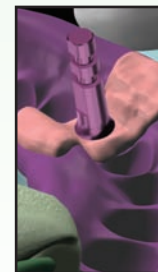


Pour the working cast in dental stone with soft tissue replication material

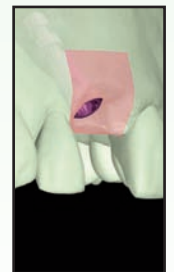


Remove the impression from the working cast exposing top of Implant Analog

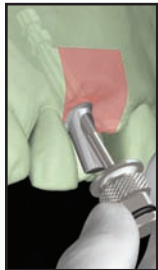

Attach the Angled Abutment to the Implant Analog in the working cast with a 1.25mmD Hex Tool

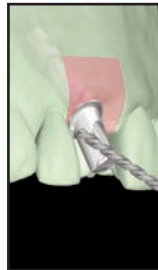

Modify angle of abutment and refine margins as required

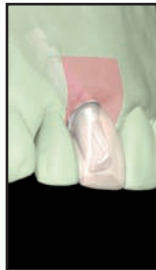

Fabricate the temporary crown



At Stage-Two Surgery, attach the prepared abutment, tighten to 30 Ncm with a Torque Wrench



Occlude screw access hole and confirm contour of restoration



Cement the provisional restoration


Common procedures for abutment placement, margin finessing, impressing and placement of restoration



Tighten to 30 Ncm with a Torque Wrench



Finesse the margins with a high-speed bur and copious irrigation prior to final impressing

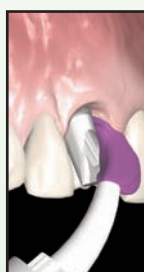


Make an impression following standard clinical procedures




Attach provisional restoration and cement onto the prepared abutment




Following standard laboratory procedures, fabricate the final restoration and cement onto the prepared abutment

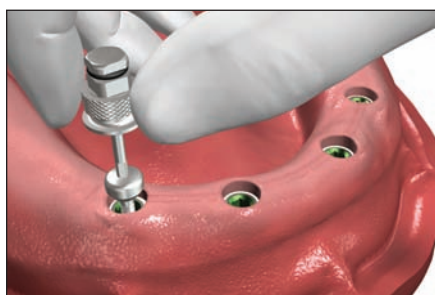


1. Unthread the healing collars from the implant using 1.25mmD Hex Tool. The top of the implants are exposed.

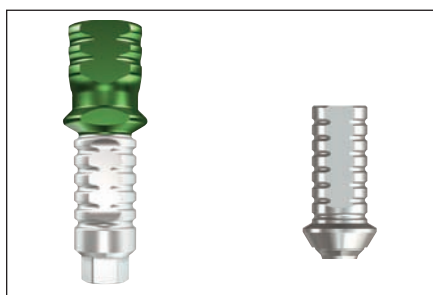

2. Use the prepared implant fixture mount as a temporary abutment (reshaping of hex at apex is essential) or select the Titanium Temporary Abutment matching the implant prosthetic platform.

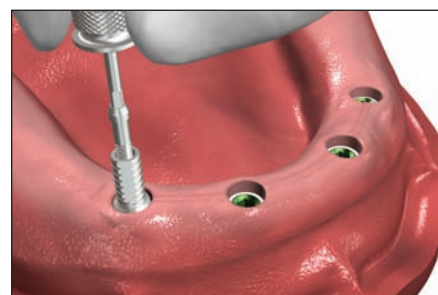

3. Place the fixation screw through the abutment and tighten in place using the appropriate 1.25mmD Hex Tool. The torque value of the fixation screw is determined by clinician. Standard torque value for fixation screws is 30Ncm.

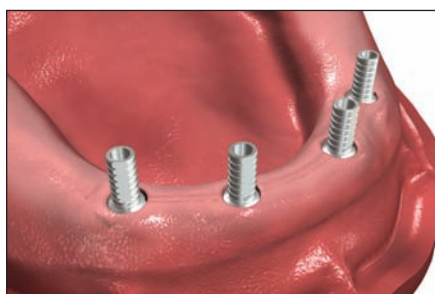


4. Confirm there is less than 40 degrees of relative divergence between abutments to be splinted together.




5. Mark position of Titanium Temporary Abutments on under surface of denture with indelible marker.




6. Use acrylic trimmer to drill holes through denture corresponding to marked positions.

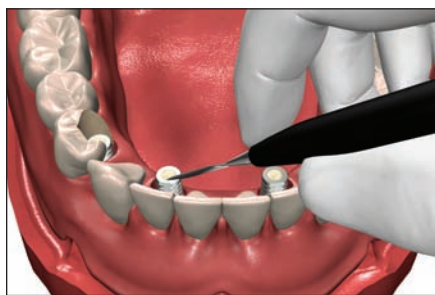


7. Seat denture over Titanium Temporary Abutment to confirm passive placement. Check occlusion with opposing arch. Occlude screw access hole with cotton pellet.

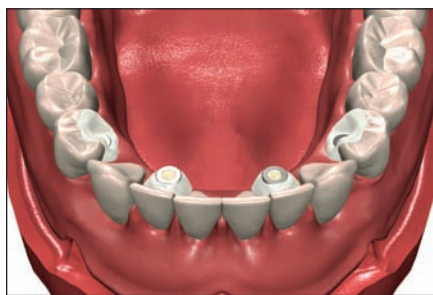

8. Use autopolymerizing acrylic or similar light cure material to attach prosthesis to Titanium Temporary Abutments. Allow to set prior to removal from mouth.

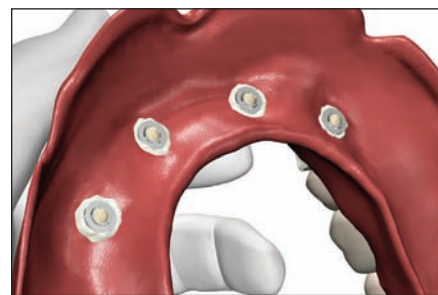


9. Remove prosthesis from the mouth, flip over to expose tissue side. Use additional autopolymerizing resin to fill in any voids around the processed abutments.

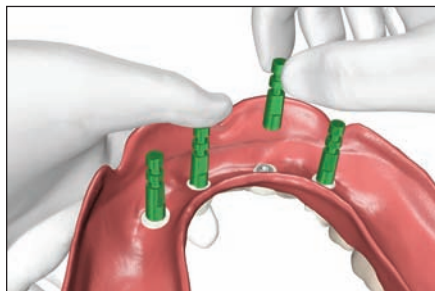

10. Optional step: Attach implant analogs to base of denture. Pour a soft tissue working cast to prevent engagement with undercuts. Articulate with opposing arch.



11. Remove denture and trim the prosthesis to allow access for patient cleaning. The prosthesis is now a fixed restoration, removable by clinician only.

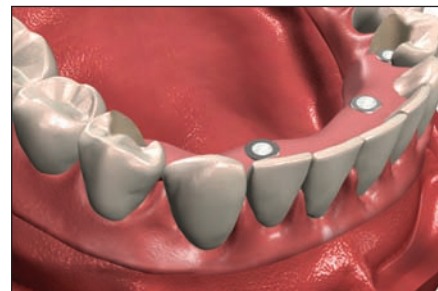


12. Attach prosthesis to implants and tighten in place using the appropriate tooling. Occlude screw access holes with material of choice.




3.0mmD Platform



3.5mmD Platform

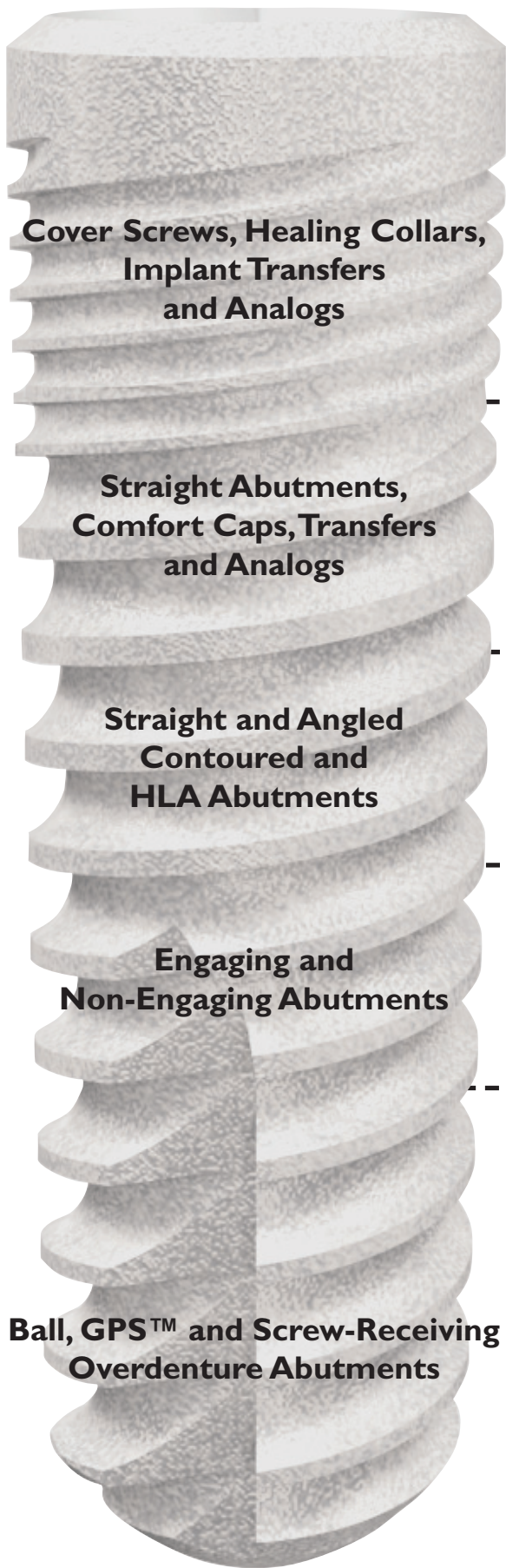


Cover Screws, Healing Collars, Implant Transfers and Analogs

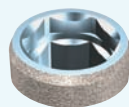

Fixture-Mount

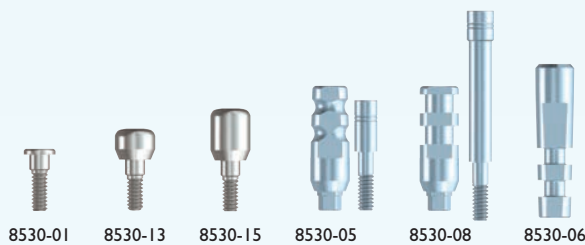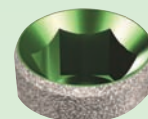


8530-01 8530-13 8530-15 8530-05 8530-08 8530-06

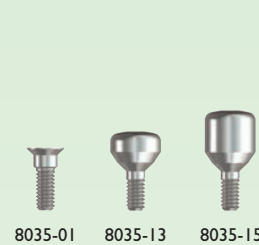


8035-01 8035-13 8035-15

Straight Abutments, Comfort Caps, Transfers and Analogs

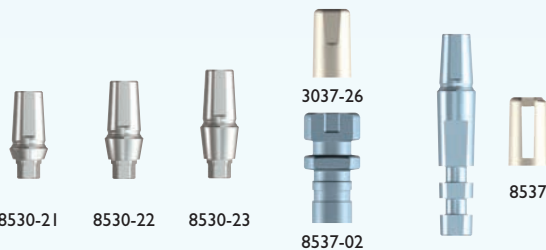

8530-21 8530-22 8530-23 3037-26 8537 8537-02

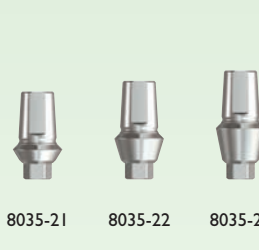

8035-21 8035-22 8035-23

Straight and Angled Contoured and HLA Abutments




8530-31 8530-32 8530-31C 8530-32C 8530-51 8530-52 8530-51-15C 8530-52-15C 8530-51-08C 8530-52-08C




8035-31 8035-32 8035-31C 8035-32C 8035-51 8035-52

Engaging and Non-Engaging Abutments

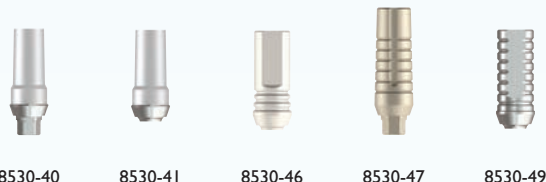

8530-40 8530-41 8530-46 8530-47 8530-49

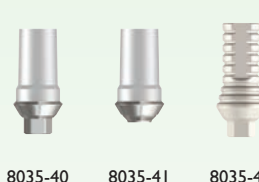

8035-40 8035-41 8035-44

Ball, GPS™ and Screw-Receiving Overdenture Abutments



8530-61 8530-62 8530-64




8035-72 8035-64-15 8035-64-30 8035-61 8035-91 8035-92 8035-93




Fixture-Mount

4.5mmD Platform

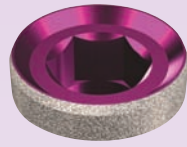


Fixture-Mount


5.7mmD Platform

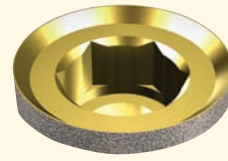


Fixture-Mount

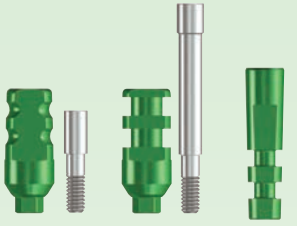

8035-05 8035-08 8035-06

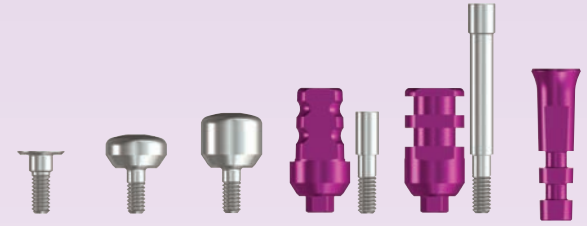


8045-01 8045-13 8045-15 8045-05 8045-08 8045-06

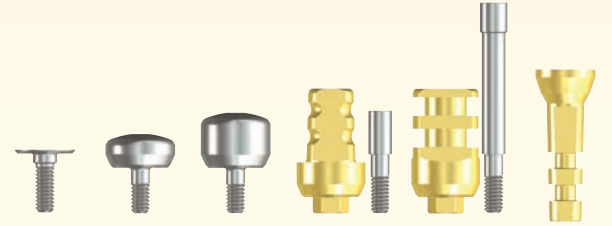


8057-01 8057-13 8057-15 8057-05 8057-08 8057-06



3047-26 3047 3047-02

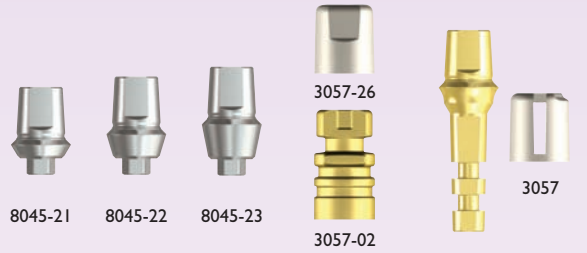

8045-21 8045-22 8045-23 3057-26 3057 3057-02

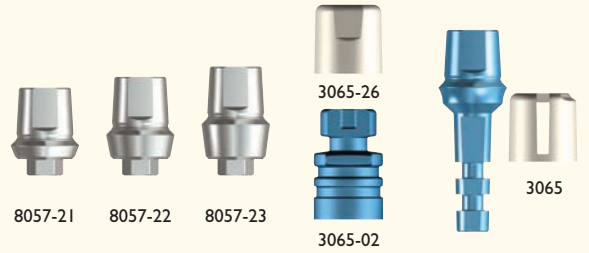

8057-21 8057-22 8057-23 3065-26 3065 3065



8035-51-15C 8035-51-08C HLA3-4 8035-52-15C 8035-52-08C HLA3-5




8045-31 8045-31C 8045-51 8045-51-15C 8045-51-08C HLA4-5 8045-32 8045-32C 8045-52 8045-52-15C 8045-52-08C HLA4-6



8057-31 8057-51 HLA5-6 8057-32 8057-52 8057-52 HLA5-7



8035-46 8035-47 8035-49

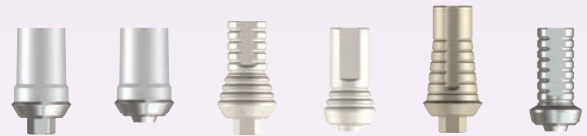

8045-40 8045-41 8045-44 8045-46 8045-47 8045-49



8057-40 8057-44 8057-46 8057-47 8057-49



8035-74



8045-72



8045-74



8057-72



8057-74



8035-62



8035-64



8045-64-15



8045-64-30



8045-61



8045-62



8045-64



8057-61



8057-62



8057-64



8035-94



8035-95



8045-91



8045-92



8045-93



8045-94



8045-95



8057-91



8057-92



8057-93



8057-94



8057-95




Cover Screws, Healing Collars, Transfers and Analogs

Surgical Cover Screws

Requires 1.25mmD Hex Tool

Code	Code [*for 6mmL Only]	Platform	Price
8530-01	-	3.0mmD	\$20.00
8035-01	8035-01S	3.5mmD	\$20.00
8045-01	8045-01S	4.5mmD	\$20.00
8057-01	8057-01S	5.7mmD	\$20.00

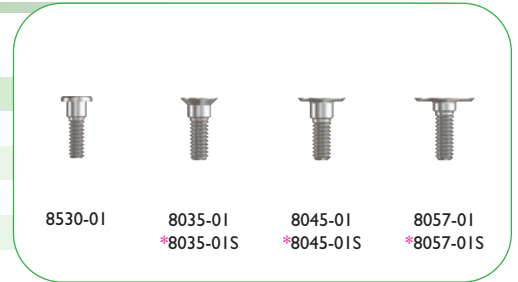


8530-01 8035-01
*8035-01S 8045-01
*8045-01S 8057-01
*8057-01S

Healing Collars

Requires 1.25mmD Hex Tool

Code	Platform	Width	Length	Price
8530-13	3.0mmD	3.8mmD	3mmL	\$30.00
8530-15	3.0mmD	3.8mmD	5mmL	\$30.00
8035-13	3.5mmD	4.8mmD	3mmL	\$30.00
8035-15	3.5mmD	4.8mmD	5mmL	\$30.00
8045-13	4.5mmD	5.8mmD	3mmL	\$30.00
8045-15	4.5mmD	5.8mmD	5mmL	\$30.00
8057-13	5.7mmD	6.6mmD	3mmL	\$30.00
8057-15	5.7mmD	6.6mmD	5mmL	\$30.00

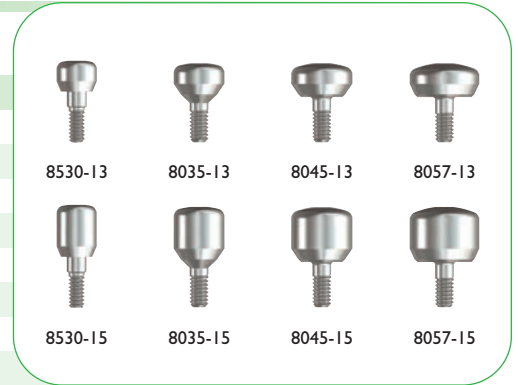


8530-13 8035-13 8045-13 8057-13
8530-15 8035-15 8045-15 8057-15

Closed-Tray Implant Transfers

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Replacement Screw	Price
8530-05	3.0mmD	1000-09	\$30.00
8035-05	3.5mmD	1000-06	\$30.00
8045-05	4.5mmD	1000-06	\$30.00
8057-05	5.7mmD	1000-06	\$30.00

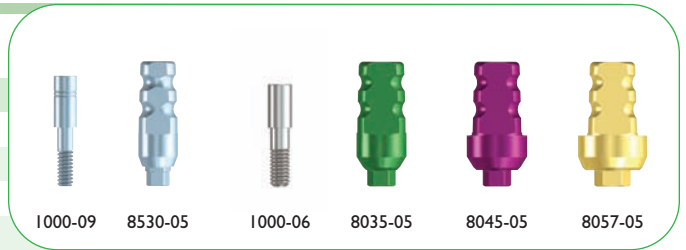


1000-09 8530-05 1000-06 8035-05 8045-05 8057-05

Open-Tray Implant Transfers

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Replacement Screw	Price
8530-08	3.0mmD	1000-26	\$30.00
8035-08	3.5mmD	1000-10	\$30.00
8045-08	4.5mmD	1000-10	\$30.00
8057-08	5.7mmD	1000-10	\$30.00

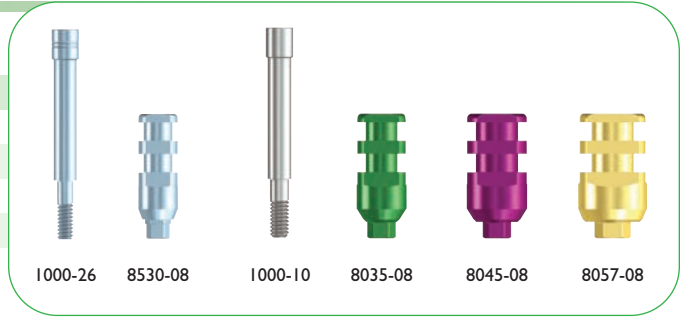

1000-26 8530-08 1000-10 8035-08 8045-08 8057-08

Implant Analogs

Code	Platform	Price
8530-06	3.0mmD	\$20.00
8035-06	3.5mmD	\$20.00
8045-06	4.5mmD	\$20.00
8057-06	5.7mmD	\$20.00




8530-06 8035-06 8045-06 8057-06

Replacement Fixation Screws

Requires 1.25mmD Hex Tool

Code	Description	Price (includes two pieces)
1000-04	For use with 3.0 hexagon platform	\$30.00
1000-09	For use with 3.0 hexagon platform	\$30.00
1000-26	For use with 3.0 hexagon platform	\$30.00
1000-01	For use with 3.5, 4.5 & 5.7mmD hexagon platforms	\$30.00
1000-06	For use with 3.5, 4.5 & 5.7mmD hexagon platforms	\$30.00
1000-10	For use with 3.5, 4.5 & 5.7mmD hexagon platforms	\$30.00

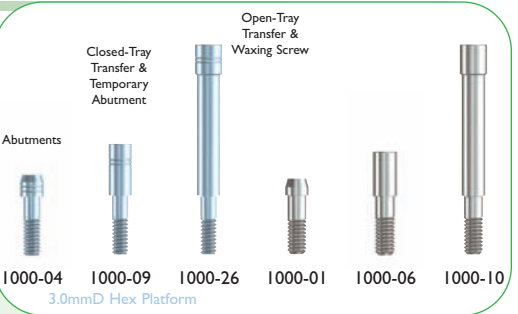

Abutments
1000-04 1000-09 1000-26 1000-01 1000-06 1000-10
3.0mmD Hex Platform
Closed-Tray Transfer & Temporary Abutment
Open-Tray Transfer & Waxing Screw




Straight Abutments, Snap-on Transfers and Lab Analogs

Straight Abutments

Includes Transfer, Comfort Cap & Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Use Analog	Price
8530-21	3.0mmD	1.0mm	8537	\$85.00
8530-22	3.0mmD	2.0mm	8537	\$85.00
8530-23	3.0mmD	3.0mm	8537	\$85.00
8035-21	3.5mmD	1.0mm	3047	\$85.00
8035-22	3.5mmD	2.0mm	3047	\$85.00
8035-23	3.5mmD	3.0mm	3047	\$85.00
8045-21	4.5mmD	1.0mm	3057	\$85.00
8045-22	4.5mmD	2.0mm	3057	\$85.00
8045-23	4.5mmD	3.0mm	3057	\$85.00
8057-21	5.7mmD	1.0mm	3065	\$85.00
8057-22	5.7mmD	2.0mm	3065	\$85.00
8057-23	5.7mmD	3.0mm	3065	\$85.00

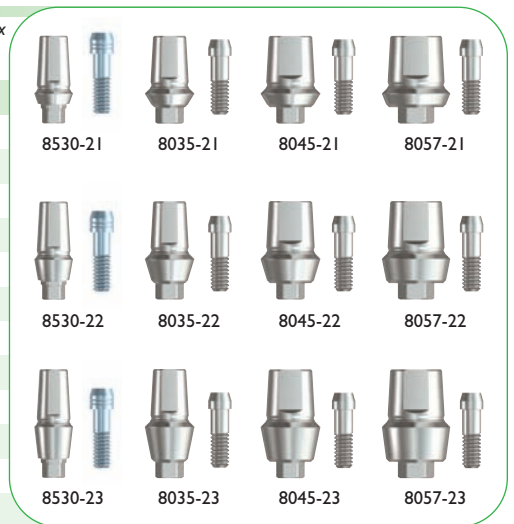


Straight Abutment Transfers

Code	Abutment Width	Price (includes two pieces)
8537-02	3.7mmD	\$30.00
3047-02	4.7mmD	\$30.00
3057-02	5.7mmD	\$30.00
3065-02	6.5mmD	\$30.00



Straight Abutment Comfort Caps

Code	Abutment Width	Price (includes two pieces)
3037-26	3.7mmD	\$30.00
3047-26	4.7mmD	\$30.00
3057-26	5.7mmD	\$30.00
3065-26	6.5mmD	\$30.00

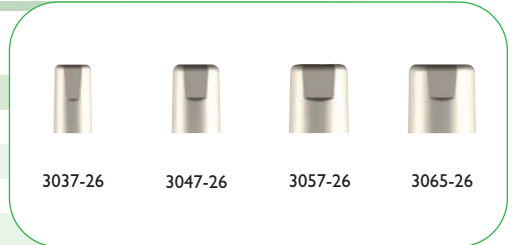

Straight Abutment Analogs

Includes Castable Coping

Code	Abutment Width	Price
8537	3.7mmD	\$35.00
3047	4.7mmD	\$35.00
3057	5.7mmD	\$35.00
3065	6.5mmD	\$35.00



Straight Abutment Castable Coping

Code	Abutment Width	Price (includes two pieces)
3037-28	3.7mmD	\$30.00
3047-28	4.7mmD	\$30.00
3057-28	5.7mmD	\$30.00
3065-28	6.5mmD	\$30.00






Straight and Angled Contoured, Titanium and Zirconia Abutments

Straight Titanium Contoured Abutments

Price includes Fixation Screw with 1.25mmD Hex


Code	Platform	Collar Heights	Replacement Screw	Price
8530-31	3.0mmD	1mm	1000-04	\$85.00
8530-32	3.0mmD	2mm	1000-04	\$85.00
8035-31	3.5mmD	1mm	1000-01	\$85.00
8035-32	3.5mmD	2mm	1000-01	\$85.00
8045-31	4.5mmD	1mm	1000-01	\$85.00
8045-32	4.5mmD	2mm	1000-01	\$85.00
8057-31	5.7mmD	1mm	1000-01	\$85.00
8057-32	5.7mmD	2mm	1000-01	\$85.00

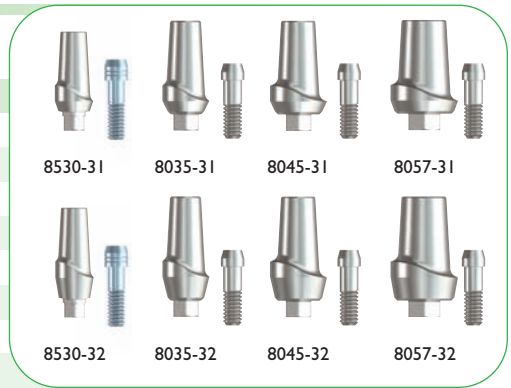

Angled Titanium Contoured Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Replacement Screw	Price
8530-51	3.0mmD	.75mm	1000-04	\$85.00
8530-52	3.0mmD	2mm	1000-04	\$85.00
8035-51	3.5mmD	.75mm	1000-01	\$85.00
8035-52	3.5mmD	2mm	1000-01	\$85.00
8045-51	4.5mmD	.75mm	1000-01	\$85.00
8045-52	4.5mmD	2mm	1000-01	\$85.00
8057-51	5.7mmD	.75mm	1000-01	\$85.00
8057-52	5.7mmD	2mm	1000-01	\$85.00

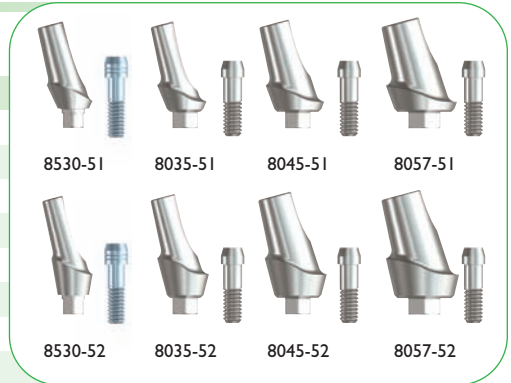

0°, 8° & 15° Zirconia Contoured Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Angulation	Collar Heights	Replacement Screw	Price
8530-31C	3.0mmD	0°	1mm	1000-04	\$120.00
8530-32C	3.0mmD	0°	2mm	1000-04	\$120.00
8530-51-08C	3.0mmD	8°	1mm	1000-04	\$120.00
8530-52-08C	3.0mmD	8°	2mm	1000-04	\$120.00
8530-51-15C	3.0mmD	15°	1mm	1000-04	\$120.00
8530-52-15C	3.0mmD	15°	2mm	1000-04	\$120.00
8530-30CB	3.0mmD	-	.5mm	1000-04	\$60.00

Code	Platform	Angulation	Collar Heights	Replacement Screw	Price
8035-31C	3.5mmD	0°	1mm	1000-01	\$120.00
8035-32C	3.5mmD	0°	2mm	1000-01	\$120.00
8035-51-08C	3.5mmD	8°	1mm	1000-01	\$120.00
8035-52-08C	3.5mmD	8°	2mm	1000-01	\$120.00
8035-51-15C	3.5mmD	15°	1mm	1000-01	\$120.00
8035-52-15C	3.5mmD	15°	2mm	1000-01	\$120.00
8035-30CB	3.5mmD	-	.5mm	1000-01	\$60.00

Code	Platform	Angulation	Collar Heights	Replacement Screw	Price
8045-31C	4.5mmD	0°	1mm	1000-01	\$120.00
8045-32C	4.5mmD	0°	2mm	1000-01	\$120.00
8045-51-08C	4.5mmD	8°	1mm	1000-01	\$120.00
8045-52-08C	4.5mmD	8°	2mm	1000-01	\$120.00
8045-51-15C	4.5mmD	15°	1mm	1000-01	\$120.00
8045-52-15C	4.5mmD	15°	2mm	1000-01	\$120.00
8045-30CB	4.5mmD	-	.5mm	1000-01	\$60.00





Engaging and Non-Engaging, Gold, Plastic and Titanium Abutments

Straight Full-Contoured (HLA) Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Replacement Screw	Price
HLA3-4	3.5mmD	0mm	1000-01	\$75.00
HLA3-5	3.5mmD	0mm	1000-01	\$75.00
HLA4-5	4.5mmD	0mm	1000-01	\$75.00
HLA4-6	4.5mmD	0mm	1000-01	\$75.00
HLA5-6	5.7mmD	0mm	1000-01	\$75.00
HLA5-7	5.7mmD	0mm	1000-01	\$75.00



Gold Engaging Cast-to Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Replacement Screw	Price
8530-40	3.0mmD	1.0mm	1000-04	\$100.00
8035-40	3.5mmD	1.0mm	1000-01	\$100.00
8045-40	4.5mmD	1.0mm	1000-01	\$100.00
8057-40	5.7mmD	1.0mm	1000-01	\$125.00



Gold Non-Engaging Cast-to Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Replacement Screw	Price
8530-41	3.0mmD	1.0mm	1000-04	\$100.00
8035-41	3.5mmD	1.0mm	1000-01	\$100.00
8045-41	4.5mmD	1.0mm	1000-01	\$100.00



Plastic Engaging Castable Abutments

Price includes Fixation Screw with 1.25mmD Hex


Code	Platform	Collar Heights	Replacement Screw	Price
8035-44	3.5mmD	0mm	1000-01	\$35.00
8045-44	4.5mmD	0mm	1000-01	\$35.00
8057-44	5.7mmD	0mm	1000-01	\$35.00



Plastic Non-Engaging Castable Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Replacement Screw	Price
8530-46	3.0mmD	0mm	1000-04	\$30.00
8035-46	3.5mmD	0mm	1000-01	\$30.00
8045-46	4.5mmD	0mm	1000-01	\$30.00
8057-46	5.7mmD	0mm	1000-01	\$30.00

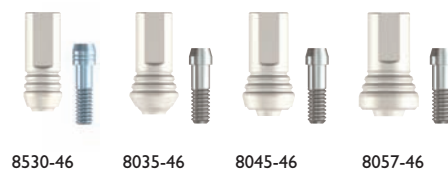

Plastic Engaging Temporary Abutments

Price includes Fixation Screw with 1.25mmD Hex


Code	Platform	Collar Heights	Replacement Screw	Price
8530-47	3.0mmD	1.0mm	1000-09	\$35.00
8035-47	3.5mmD	1.0mm	1000-06	\$35.00
8045-47	4.5mmD	1.0mm	1000-06	\$35.00
8057-47	5.7mmD	1.0mm	1000-06	\$35.00



Titanium Non-Engaging Temporary Abutments

Price includes Fixation Screw with 1.25mmD Hex

Code	Platform	Collar Heights	Replacement Screw	Price
8530-49	3.0mmD	1.0mm	1000-04	\$35.00
8035-49	3.5mmD	1.0mm	1000-01	\$35.00
8045-49	4.5mmD	1.0mm	1000-01	\$35.00
8057-49	5.7mmD	1.0mm	1000-01	\$35.00

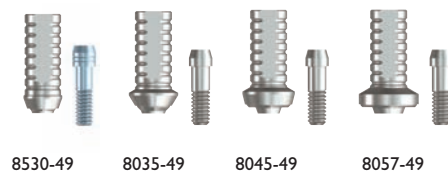



Selection Chart: Overdenture and Screw-Retained Restorations

	3.0mmD Platform	3.5mmD Platform	4.5mmD Platform	5.7mmD Platform					
Screw-Receiving Overdenture Abutments (*Includes Transfer and Comfort Cap)									
Overdenture Abutments	 8530-61 8530-62 8530-64	 8035-61 8035-62 8035-64 8035-64-15 8035-64-30	 8045-61 8045-62 8045-64 8045-64-15 8045-64-30	 8057-61 8057-62 8057-64					
Screw-Receiving Overdenture Abutment Transfer, Analog & Comfort Cap									
OD Abutment Comfort Cap Transfer and Analogs	 *4000-13	 *1000-68	 1000-61	 1000-62 1000-63					
Overdenture Screw-Retained Copings									
Overdenture Abutment Components	 1000-64 Gold	 1000-65 Titanium	 1000-66 Plastic	 1000-67 Gold/Plastic	 1000-69 Titanium				
Ball Abutments (**Includes Transfer and Comfort Cap)									
Ball Abutments		 8035-72 8035-74	 8045-72 8045-74	 8057-72 8057-74					
Ball Attachments									
Ball Abutment Components	 **1000-87	 **1000-84	 1000-85	 1000-81	 1000-82	 1000-88	 1000-83	 1000-86	
GPS Abutments (***)Includes Transfer, Abutment Housing Kit and Comfort Cap)									
GPS Abutments	 8530-91 8530-92 8530-93 8530-94	 8035-91 8035-92 8035-93 8035-94 8035-95	 8045-91 8045-92 8045-93 8045-94 8045-95	 8057-91 8057-92 8057-93 8057-94 8057-95					
GPS Attachments									
GPS Abutment Components	 ***4500-26	 ***1000-94	 1000-95	 1000-92R	 1000-92G	 1000-92E	 ***1000-91 1000-93	 1000-98S	 1000-96

Abutments are compatible with Zimmer, BioHorizons (internal) and MIS Implants with 3.5mm, 4.5mmD and 5.7mmD prosthetic platforms




Screw-Receiving, Ball and GPS Abutments

Screw-Receiving Overdenture Abutments

Includes Transfer & Comfort Cap - requires 1.25mmD Hex Tool

Platform	Collar Heights			Price
	1mmL	2mmL	4mmL	
3.0mmD	-	8530-62	8530-64	\$85.00
3.5mmD	8035-61	8035-62	8035-64	\$85.00
4.5mmD	8045-61	8045-62	8045-64	\$85.00
5.7mmD	8057-61	8057-62	8057-64	\$85.00

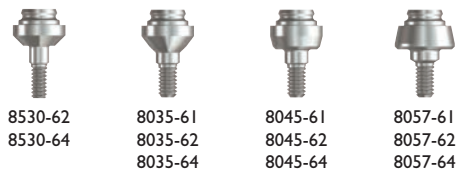

Ball Abutments

Includes Transfer & Comfort Cap - requires 1.25mmD Hex Tool


Platform	Collar Heights		Price
	2mmL	4mmL	
3.5mmD	8035-72	8035-74	\$85.00
4.5mmD	8045-72	8045-74	\$85.00
5.7mmD	8057-72	8057-74	\$85.00



GPS Abutments

Includes Transfer, Comfort Cap, Metal Housing, Clear Male Liner & Spacer - requires 1.25mmD Hex Tool


Platform	Collar Heights					Price
	1mmL	2mmL	3mmL	4mmL	5mmL	
3.0mmD	8530-91	8530-92	8530-93	8530-94	-	\$100.00
3.5mmD	8035-91	8035-92	8035-93	8035-94	8035-95	\$100.00
4.5mmD	8045-91	8045-92	8045-93	8045-94	8045-95	\$100.00
5.7mmD	8057-91	8057-92	8057-93	8057-94	8057-95	\$100.00

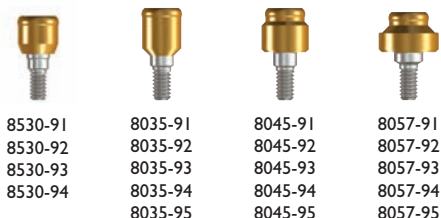

Screw-Receiving Angled Overdenture Abutments

Includes Transfer & Comfort Cap - requires 1.25mmD & 2.5mmD Hex Tool


Code	Platform	Description	Price
8035-64-15	3.5mmD	15° Angled Abutment	\$100.00
8045-64-15	4.5mmD	15° Angled Abutment	\$100.00
8035-64-30	3.5mmD	30° Angled Abutment	\$100.00
8045-64-30	4.5mmD	30° Angled Abutment	\$100.00

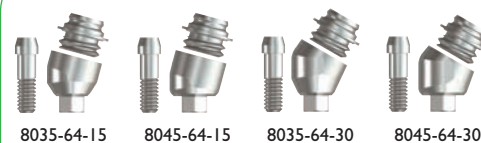

Screw-Receiving Angled GPS Abutments

Includes Transfer, Comfort Cap, Metal Housing, Clear Male Liner & Spacer - requires 1.25mmD Hex Tool


Code	Platform	Description	Price
8035-94-15	3.5mmD	15° Angled Abutment	\$120.00
8045-94-15	4.5mmD	15° Angled Abutment	\$120.00
8035-94-30	3.5mmD	30° Angled Abutment	\$120.00
8045-94-30	4.5mmD	30° Angled Abutment	\$120.00

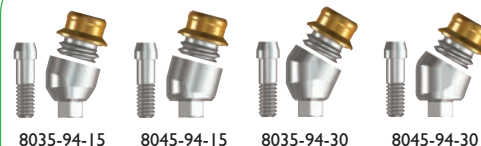

Replacement Tops for Angled Abutments

Includes Transfer & Comfort Cap

Code	Description	Price
1000-60	Replacement Screw Top for Angled Base <small>[Requires 2.5mmD Hex Tool]</small>	\$35.00
1000-70	Replacement Ball Top for Angled Base <small>[Requires 1.25mmD Hex Tool]</small>	\$35.00
1000-90	Replacement GPS Top for Angled Base <small>[Requires 1.25mmD Hex Tool]</small>	\$35.00

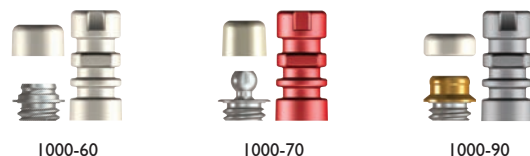




Ancillary Components for Overdenture and Screw-Retained Restorations


Screw-Receiving Overdenture Abut. Comfort Cap, Transfers & Analog Screws require 1.25mmD Hex Tool

Code	Prosthetic Platform	Description	Price
1000-61	5.0mmD	Aluminum Analog	\$20.00
1000-62	5.0mmD	Open-Tray Transfer plus Fixation Screw	\$30.00
1000-63	5.0mmD	Closed-Tray Transfer plus Fixation Screw	\$30.00
1000-68	5.0mmD	Plastic Abutment Transfer - 2ea	\$30.00
4000-13	5.0mmD	Plastic Comfort Cap - 2ea	\$30.00
1000-14	-	Replacement Open-Tray Fixation Screw - 2ea	\$30.00
1000-15	-	Replacement Closed-Tray Fixation Screw - 2ea	\$30.00

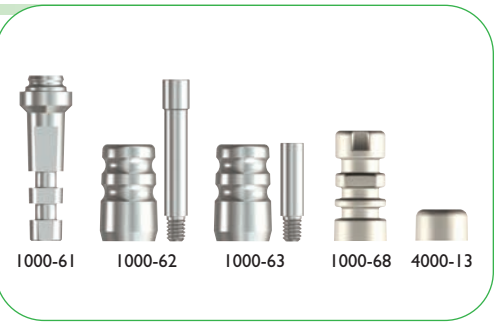


Overdenture Screw-Retained Abutment Components Includes Fixation Screw [1000-05] with 1.25mmD Hex

Code	Prosthetic Platform	Description	Price
1000-05	-	Coping Replacement Fixation Screw - 2ea	\$30.00
1000-64	5.0mmD	Gold Coping	\$80.00
1000-65	5.0mmD	Titanium Coping	\$40.00
1000-66	5.0mmD	Plastic Coping	\$35.00
1000-67	5.0mmD	Gold & Plastic Coping	\$80.00
1000-69	5.0mmD	Titanium Temporary Coping	\$35.00

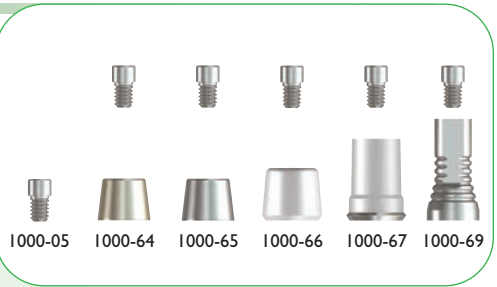

Ball Abutment Components

Code	Description	Price
1000-81	Ball Abutment Housing - 2ea	\$30.00
1000-82	Nylon Liner - 2ea	\$15.00
1000-83	Nylon Liner Transfer - 2ea	\$15.00
1000-84	Ball Abutment Transfer - 2ea	\$30.00
1000-85	Ball Abutment Analog	\$20.00
1000-86	Ball Abutment Tools	\$65.00
1000-87	Ball Abutment Comfort Cap - 2ea	\$30.00
1000-88	Ball Abutment Housing & Nylon Liner	\$22.00

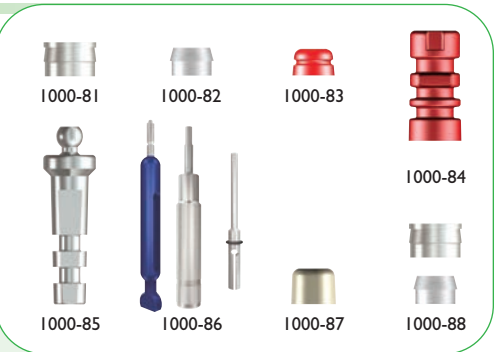

GPS Abutment Transfer, Analogs and Alignment Tools

Code	Description	Price
1000-94	Impression Coping - 2ea	\$30.00
1000-95	Abutment Analog 4mmD - 2ea	\$32.00
1000-97	Abutment Parallel Post - 4ea	\$16.00
1000-10	Guide Screw for Legacy - 2ea	\$30.00
1000-26	Guide Screw for Legacy 3.0mmD Platform - 2ea	\$30.00
4500-26	Comfort Cap - 2ea	\$10.00
1000-99	Plastic Abutment Angle Measure Guide	\$5.00



GPS Abutment Components

Code	Retention Force	Description	Price
1000-92C	4.5lbs	Replacement Male Clear - 4ea	\$16.00
1000-92G	3.0lbs	Replacement Male Green - 4ea	\$16.00
1000-92E	2.0lbs	Replacement Male Orange - 4ea	\$16.00
1000-92R	1.0lb	Replacement Male Red - 4ea	\$16.00
1000-91	-	GPS Abutment Housing Kit - 1ea	\$10.00
1000-93	-	Black Processing Male - 4ea	\$16.00
1000-98S	Assorted	Abutment Processing Cap/Male Kit - 1ea	\$20.00



GPS Instrumentation


Code	Description	Price
1000-96	Core Tool - Male Removal & Seating Tool	\$75.00





Tapered Hex Tools


Code	Description	Price
THT1.25	1.25mmD Tapered Hex Tool Long	\$40.00
THT1.25S	1.25mmD Tapered Hex Tool Short	\$40.00

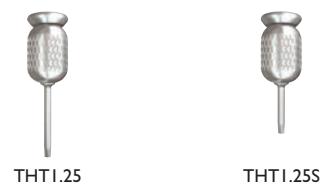

Straight Hex Tools

Square connection for use with Ratchet or Torque Wrench

Code	Description	Price
HT1.25	1.25mmD Straight Hex Tool Long (for use with Torque Wrench)	\$40.00
HT1.25S	1.25mmD Straight Hex Tool Short (for use with Torque Wrench)	\$40.00
HT2.0	2.0mmD Hex Tool Long	\$40.00
HT2.0S	2.0mmD Hex Tool Short	\$40.00
HT2.5	2.5mmD Hex Tool Long (with retentive O-ring)	\$40.00
HT2.5S	2.5mmD Hex Tool Short (with retentive O-ring)	\$40.00
HT3.0	3.0mmD Hex Tool Long	\$40.00
HT3.0S	3.0mmD Hex Tool Short	\$40.00

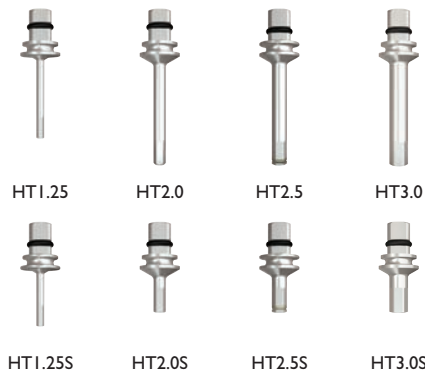

Hex Drills


Latch head connection for use with Surgical Motor

Code	Description	Price
HD1.25	1.25mmD Hex Drill	\$40.00
HD2.0	2.0mmD Hex Drill	\$40.00
HD2.5	2.5mmD Hex Drill (with retentive O-ring)	\$40.00
HD3.0	3.0mmD Hex Drill	\$40.00



Ancillary Instruments

Code	Description	Price
R	Ratchet	\$125.00
TW30	30Ncm Torque Wrench	\$125.00
TW35	35Ncm Torque Wrench	\$125.00
SSH2.5	2.5mmD Stainless Steel Handle (with retentive O-ring)	\$100.00
FMD	Fixture Mount Drill	\$80.00
PAR	Paralleling Tool 2.8/2.3	\$15.00
LDRILL	Locator Drill	\$30.00
DE	Drill Extender	\$40.00
TWM	Multi-Setting Torque Wrench (15, 20, 25, 30, 35, 40, 45, 50 & 60Ncm)	\$395.00
ABH	Analog Holder	\$85.00

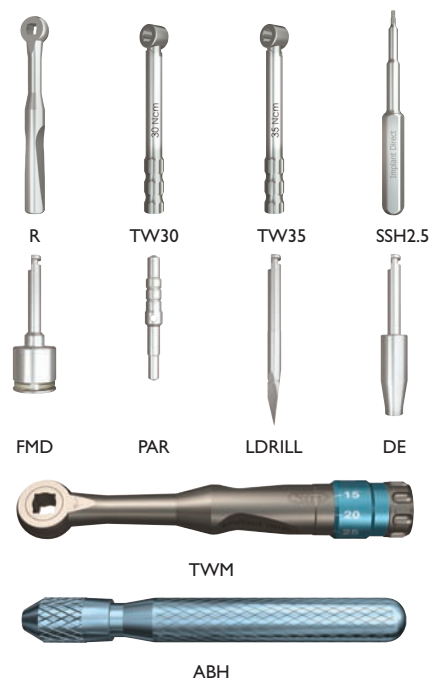

Surgical Trays

Code	Description	Price
ST	Empty Surgical Kit for Legacy System Implants	\$250.00
BST	Basic Surgical Kit for Legacy System Implants	\$1265.00
CST	Complete Surgical Kit for Legacy System Implants	\$2565.00





3.2mmD Implants

Code	Description	Price
SD2.3	2.3/2.0mmD Step Drill with 6-16 markings	\$60.00
SD2.3S	2.3/2.0mmD Step Drill with 6-10 markings	\$60.00
SD2.8	2.8/2.3mmD Step Drill with 6-16 markings	\$60.00
SD2.8S	2.8/2.3mmD Step Drill with 6-10 markings	\$60.00
CBD3.2	Cortical Bone Drill for 3.0 and 3.2mmD Implants	\$60.00



3.7mmD Implants

Code	Description	Price
SD3.4	3.4/2.8mmD Step Drill with 6-16 markings	\$60.00
SD3.4S	3.4/2.8mmD Step Drill with 6-10 markings	\$60.00
CBD3.7	Cortical Bone Drill for 3.7mmD Implants	\$60.00



4.2mmD Implants

Code	Description	Price
SD3.8	3.8/3.4mmD Step Drill with 6-16 markings	\$60.00
SD3.8S	3.8/3.4mmD Step Drill with 6-10 markings	\$60.00
CBD4.2	Cortical Bone Drill for 4.2mmD Implants	\$60.00



4.7mmD Implants

Code	Description	Price
SD4.4	4.4/3.8mmD Step Drill with 6-16 markings	\$60.00
SD4.4S	4.4/3.8mmD Step Drill with 6-10 markings	\$60.00
CBD4.7	Cortical Bone Drill for 4.7mmD Implants	\$60.00




5.2mmD Implants

Code	Description	Price
SD4.8	4.8/4.4mmD Step Drill with 6-16 markings	\$60.00
SD4.8S	4.8/4.4mmD Step Drill with 6-10 markings	\$60.00
CBD5.2	Cortical Bone Drill for 5.2mmD Implants	\$60.00



5.7 and 7.0mmD Implants

Code	Description	Price
D5.1	5.1mmD Straight Drill with 6-16 markings	\$60.00
D5.1S	5.1mmD Straight Drill with 6-10 markings	\$60.00
SD5.4	5.4/4.8mmD Step Drill with 6-16 markings	\$60.00
SD5.4S	5.4/4.8mmD Step Drill with 6-10 markings	\$60.00
CBD5.7	Cortical Bone Drill for 5.7mmD Implants	\$60.00
CBD7.0	Cortical Bone Drill for 7.0mmD Implants	\$60.00

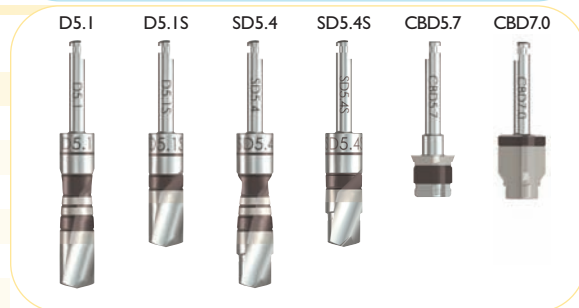



Drill Sequence

3.2mmD Straight /Tapered Implant

3.0/3.2



Soft Bone Final



SD2.3
2.3/2.0mmD

Dense Bone Final



SD2.8
2.8/2.3mmD

Optional
Bone Contour



CBD3.2

3.7mmD Straight /Tapered Implant

3.5/3.7



SD2.3
2.3/2.0mmD

Soft Bone Final



SD2.8
2.8/2.3mmD

Dense Bone Final



SD3.4
3.4/2.8mmD

Optional
Bone Contour



CBD3.7

4.2mmD Straight /Tapered Implant

3.5/4.2



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD

Soft Bone Final



SD3.4
3.4/2.8mmD

Dense Bone Final



SD3.8
3.8/3.4mmD

Optional
Bone Contour



CBD4.2

4.7mmD Straight /Tapered Implant

4.5/4.7



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD



SD3.4
3.4/2.8mmD

Soft Bone Final



SD3.8
3.8/3.4mmD

Dense Bone Final



SD4.4
4.4/3.8mmD

Optional
Bone Contour



CBD4.7

5.2mmD Straight /Tapered Implant

4.5/5.2



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD



SD3.4
3.4/2.8mmD

Optional



SD3.8
3.8/3.4mmD

Soft Bone Final



SD4.4
4.4/3.8mmD

Dense Bone Final



SD4.8
4.8/4.4mmD

Optional
Bone Contour



CBD5.2

5.7mmD & 7.0mmD Straight /Tapered Implant

5.7/5.7

5.7/7.0



SD2.3
2.3/2.0mmD

Optional



SD2.8
2.8/2.3mmD



SD3.4
3.4/2.8mmD

Optional



SD3.8
3.8/3.4mmD



SD4.4
4.4/3.8mmD

Optional



SD4.8
4.8/4.4mmD

Soft Bone Final



D5.1
5.1mmD

Dense Bone Final



SD5.4
5.4/4.8mmD

Optional
Bone Contour
for 5.7mmD



CBD5.7


Optional
Bone Contour
for 7.0mmD



CBD7.0



Complete Surgical Tray Contents

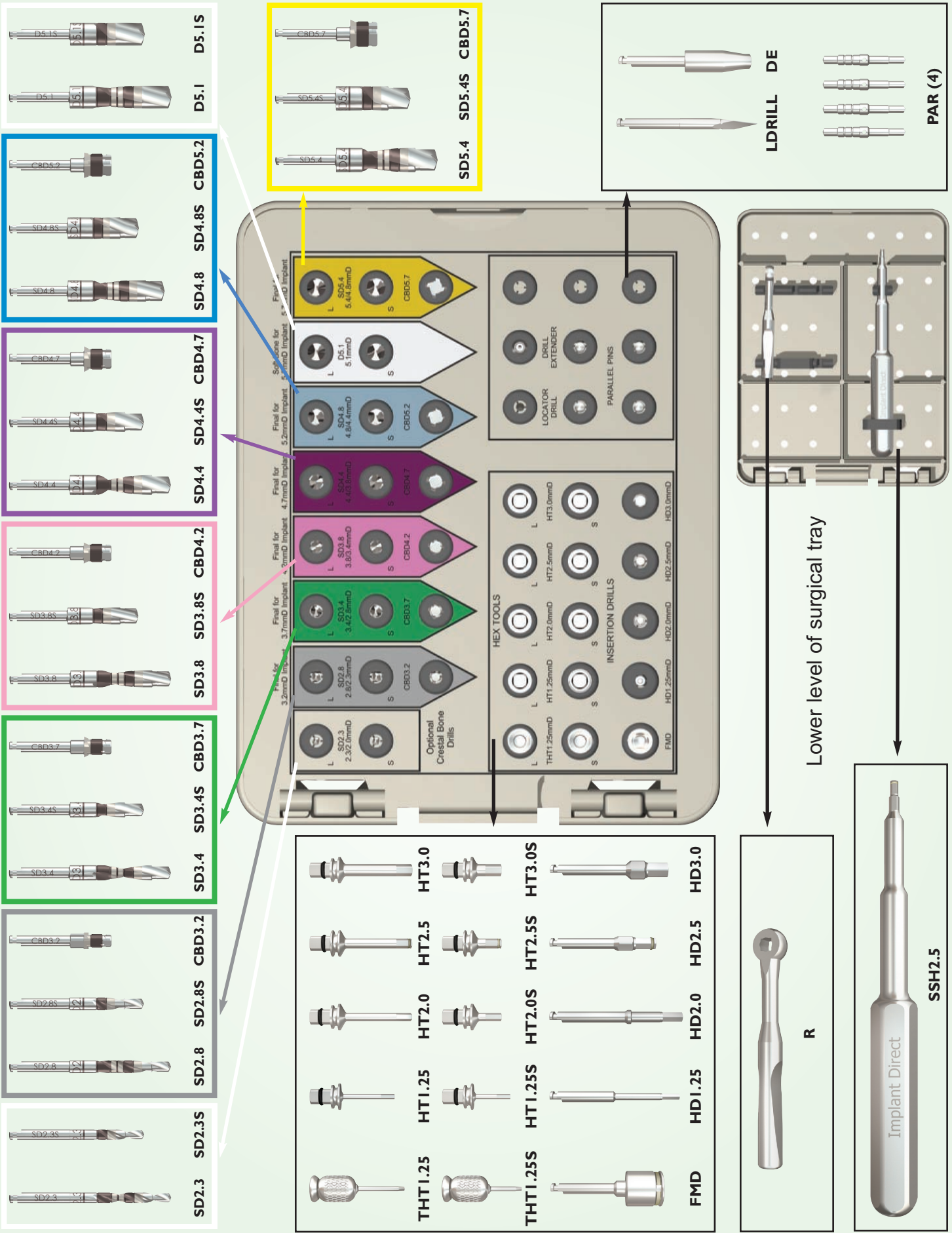



Lid



Inside



Empty Tray

Surgical Tray Contents			Complete Kit	Basic Kit	Conversion Kit
Code	Description	Price	\$2,565.00	\$1,265.00	\$930.00
ST	Empty Surgical Tray	\$250.00	X	X	X
SD2.3	2.3/2.0mmD Step Drill-Long	\$60.00	X	X	
SD2.3S	2.3/2.0mmD Step Drill-Short	\$60.00	X		
SD2.8	2.8/2.3mmD Step Drill-Long	\$60.00	X	X	
SD2.8S	2.8/2.3mmD Step Drill-Short	\$60.00	X		
CBD3.2	3.2mmD Cortical Drill	\$60.00	X		X
SD3.4	3.4/2.8mmD Step Drill-Long	\$60.00	X	X	
SD3.4S	3.4/2.8mmD Step Drill-Short	\$60.00	X		
CBD3.7	3.7mmD Cortical Drill	\$60.00	X		X
SD3.8	3.8/3.4mmD Step Drill-Long	\$60.00	X	X	X
SD3.8S	3.8/3.4mmD Step Drill-Short	\$60.00	X		X
CBD4.2	4.2mmD Cortical Drill	\$60.00	X		X
SD4.4	4.4/3.8mmD Step Drill-Long	\$60.00	X	X	
SD4.4S	4.4/3.8mmD Step Drill-Short	\$60.00	X		
CBD4.7	4.7mmD Cortical Drill	\$60.00	X		X
SD4.8	4.8/4.4mmD Step Drill-Long	\$60.00	X	X	X
SD4.8S	4.8/4.4mmD Step Drill-Short	\$60.00	X		X
CBD5.2	5.2mmD Cortical Drill	\$60.00	X		X
D5.1	5.1mmD Straight Drill-Long	\$60.00	X	X	
D5.1S	5.1mmD Straight Drill-Short	\$60.00	X		
SD5.4	5.4/4.8mmD Step Drill-Long	\$60.00	X	X	
SD5.4S	5.4/4.8mmD Step Drill-Short	\$60.00	X		
CBD5.7	5.7mmD Cortical Drill	\$60.00	X		X
Hex Tools					
THT1.25	1.25mmD Tapered Hex Tool-Long	\$40.00	X	X	
THT1.25S	1.25mmD Tapered Hex Tool-Short	\$40.00	X		
HT1.25	1.25mmD Straight Hex Tool-Long	\$40.00	X	X	
HT1.25S	1.25mmD Straight Hex Tool-Short	\$40.00	X		
HT2.0	2.0mmD Straight Hex Tool-Long	\$40.00	X	X	
HT2.0S	2.0mmD Straight Hex Tool-Short	\$40.00	X		
HT2.5	2.5mmD Straight Hex Tool-Long	\$40.00	X	X	
HT2.5S	2.5mmD Straight Hex Tool-Short	\$40.00	X		
HT3.0	3.0mmD Straight Hex Tool-Long	\$40.00	X	X	
HT3.0S	3.0mmD Straight Hex Tool-Short	\$40.00	X		
HD1.25	1.25mmD Hex Drill	\$40.00	X	X	
HD2.0	2.0mmD Hex Drill	\$40.00	X		X
HD2.5	2.5mmD Hex Drill	\$40.00	X	X	
HD3.0	3.0mmD Hex Drill	\$40.00	X		X
Ancillary Instrumentation					
FMD	Fixture Mount Drill	\$80.00	X		
LDRILL	Locator Drill	\$30.00	X	X	
DE	Drill Extender	\$40.00	X	X	
R	Ratchet	\$125.00	X	X	
PAR (x4)	Paralleling Tool 2.8/2.3	\$60.00	X	X	
SSH2.5	2.5mmD Stainless Steel Handle	\$100.00	X		



www.implantdirect.com
 Implant Direct Sybron International
 1-818-444-3333 • 1-888-NIZNICK



"Order Online for High Quality Products at Factory-Direct Prices"

Gerald Niznick, DMD,MSD, Founder and President of Implant Direct

Legacy™3 All-in-1 Packaging

INCLUDES: Cover Screw, Healing Collar, Transfer & Straight Preparable Abutment. Surgically and Prosthetically Compatible with Zimmer's Tapered Screw-Vent®

Legacy3 USA price = \$175

*Tapered Screw-Vent USA price = \$584

* Includes similar components provided with Legacy3 Implant



Three Implant Designs and Packaging Options
Original Conical Connection (GNiznick US Pat.#4,960,381)
 Implants include Cover Screw and Healing Collar


Select Implant based on Price, Packaging and Thread Design

Legacy1: \$125 Original "V" Thread with Plastic Carrier

Legacy2: \$150 includes Temporary Abutment/Transfer

Legacy3: \$175 includes Preparable Abutment/Transfer


Legacy1
 Standard "V" Threads
 Matches Screw-Vent

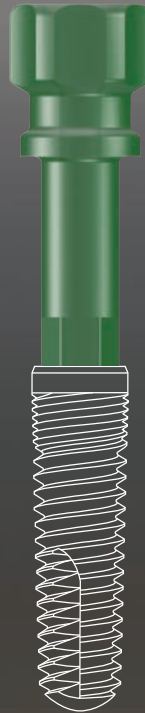

Legacy2
 Spiral Threads
 for Increased Stability



Legacy3
 Buttress Threads
 for Increased Surface

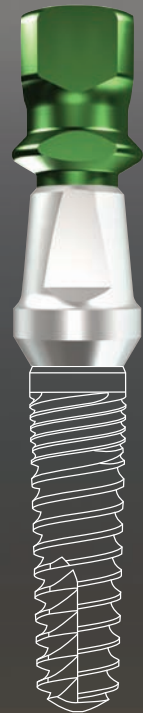

Prosthetic compatibility with Screw-Vent, BioHorizons® and MIS® Implants

INTRODUCTORY LEGACY OFFER

Try a 13mmL Legacy1, Legacy2 and Legacy3 Implant at no charge with purchase of required instruments

US List Price for 3 Implants Abutments and Prosthetic Components

Legacy = \$620

Zimmer Dental = \$1797



3.2mm Legacy3, 3.7mm Legacy1, 4.2mm Legacy2, 1 Straight Snappy and 1 Contoured Abutment with Fixation Screws



Compare US List Price	Straight Contoured	HLA Abutment	15° Angled Contoured	Straight Snap-On	Gold/Plastic	Zirconia Abutment	Plastic Temporary Abutment	Ball Attachment	GPS™ Attachment	Multiple-Unit w / Cap & Transfer	Straight	Angled
Implant Direct	\$85	\$75	\$85	\$85	\$100	\$120	\$35	\$107	\$100	\$85	\$100	\$100
Zimmer Dental	\$158	\$152	\$179	\$203	\$180	\$232	\$52	\$168	\$166	\$167	\$167	N/A
Cost Savings	46%	51%	53%	58%	44%	48%	33%	36%	40%	49%	49%	N/A

FREE Introductory Offers at
www.implantdirect.com

Terms and Conditions Apply

List prices as of 4/23/10

Screw-Vent® is a registered Trademark of Zimmer Dental

10K-12-10-0041