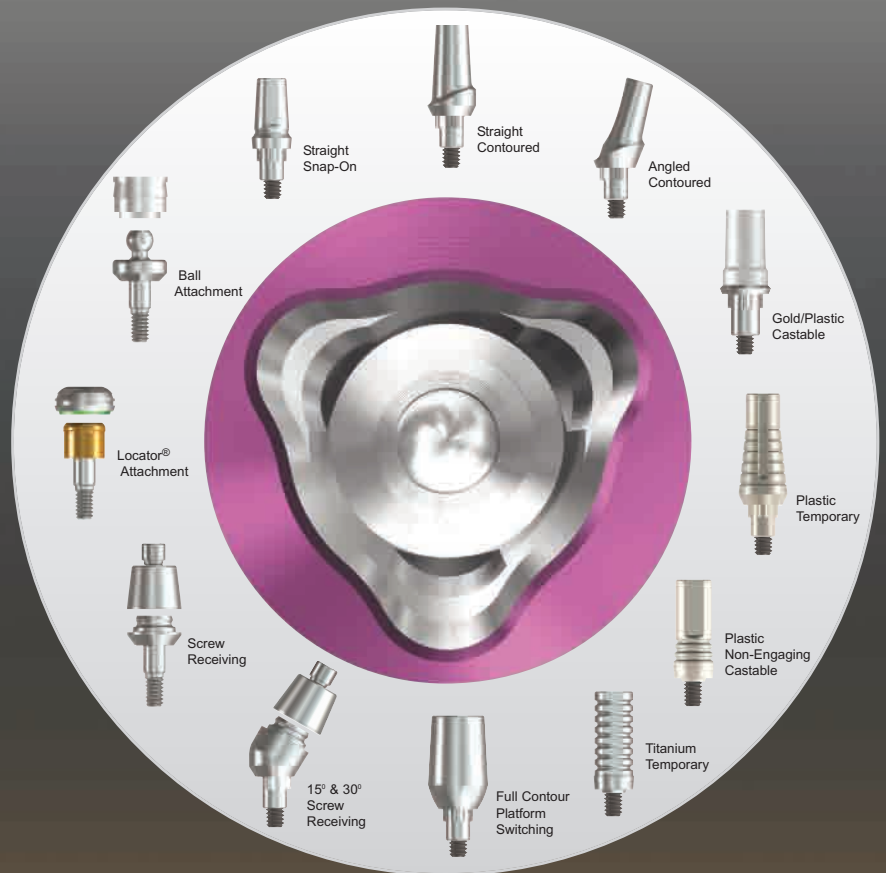




Implant
Direct™

simply smarter.

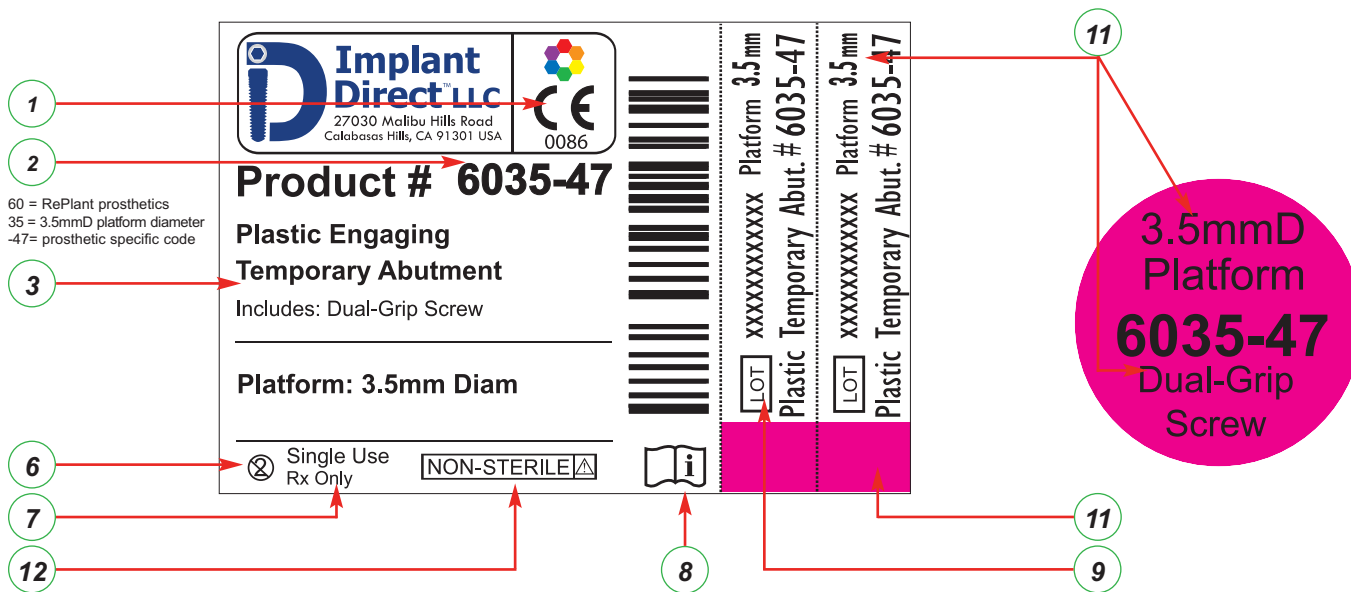
RePlant® Prosthetic Catalog



ReActive™



RePlant *Prosthetic* Vial and Cap Label



- | | |
|--|---|
| | 1) Implant direct products carry the CE mark and fulfill the requirements of the Medical Device Directive 93/42/EEC |
| | 2) Product code or reference number |
| | 3) Product description |
| | 4) Use before expiration date (XXXXXXX) |
| | 5) Gamma sterilization symbol |
| | 6) Single-use only symbol |
| | 7) Caution: US law restricts the sale of these devices by or on the order of a dentist or physician |
| | 8) Refer to Instructions For Use (IFU) |
| | 9) Lot number (XXXXXXXXXX) |
| | 10) Implant body diameter color reference on surgical tray |
| | 11) Implant prosthetic platform diameter (written and color coded) |
| | 12) Non-sterile |



The RePlant, RePlus and ReActive Implants are mounted on color-coded fixture-mounts, designed for use as implant level transfers or shortened to function as an abutment. A cover screw is included in each vial. The RePlus fixture-mount is provided with a snap-on comfort cap.



Prepare surgical site with double cutting step drill



All of Implant Direct's implants come in a sterile vial



Each label is color-coded with two tear-off labels



Each cap is color coded by platform diameter of implant



Drop implant assembly out of vial onto a sterile surface



Remove implant assembly with a fixture-mount drill or ratchet



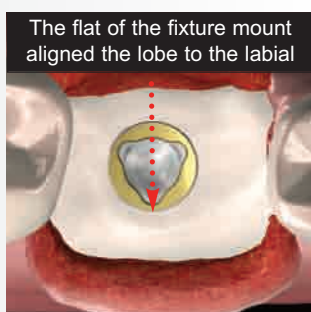
Insert implant into osteotomy



Index fixture-mount flat to labial



Remove fixture-mount




The flat of the fixture mount aligned the lobe to the labial




Remove the holder containing the cover screw & comfort cap



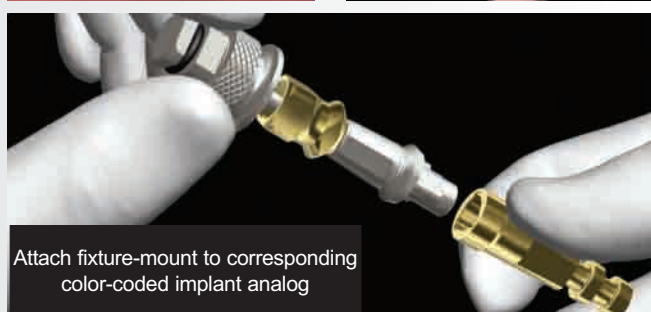Snap-on comfort cap is located in the top portion of holder



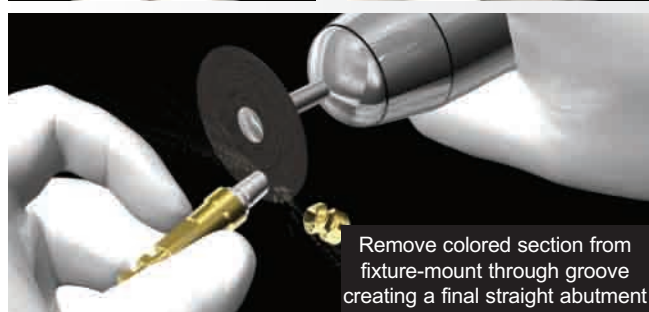Attach fixture-mount to corresponding color-coded implant analog




Remove colored section from fixture-mount through groove creating a final straight abutment



Use 30Ncm in hard bone. In soft bone, only use finger pressure



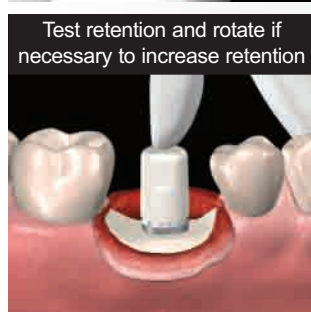Check for occlusal clearance prior to attaching comfort cap



Test retention and rotate if necessary to increase retention



Suture soft tissue around comfort cap



Selection Chart: Implant and Straight Abutments, Transfers and Analogs

3.5mmD Platform		4.3mmD Platform		5.0mmD Platform		6.0mmD Platform	
 Carrier/Transfer Implant Transfer 6035-05 6035-08		 Carrier/Transfer Implant Transfer 6043-05 6043-08		 Carrier/Transfer Implant Transfer 6050-05 6050-08		 Carrier/Transfer Implant Transfer 6060-05 6060-08	
 Implant Analog 6035-06		 Implant Analog 6043-06		 Implant Analog 6050-06		 Implant Analog 6060-06	
 Straight Abutment		 RePlus Fixture Mount		 Straight Abutment		 RePlus Fixture Mount	
 Straight Abutment		 Straight Abutment		 Straight Abutment		 Straight Abutment	
 Straight Abutment Transfer 3037-02		 Straight Abutment Transfer 3047-02		 Straight Abutment Transfer 3057-02		 Straight Abutment Transfer 3065-02	
 Comfort Cap 3037-26		 Castable Coping 3037-28		 Comfort Cap 3047-26		 Castable Coping 3047-28	
 Comfort Cap 3057-26		 Castable Coping 3057-28		 Comfort Cap 3065-26		 Castable Coping 3065-28	
 Straight Abutment Analog 3037 incl. Castable Coping		 Straight Abutment Analog 3047 incl. Castable Coping		 Straight Abutment Analog 3057 incl. Castable Coping		 Straight Abutment Analog 3065 incl. Castable Coping	




Straight Abutments, Snap-on Transfers and Lab Analogs

Straight Abutments

Includes Transfer, Comfort Cap & Fixation Screw with Dual-Grip™

Code	Collar Heights	Platform	Price
6035-21	1mm	3.5mmD	\$85.00
6035-22	2mm	3.5mmD	\$85.00
6035-23	3mm	3.5mmD	\$85.00
6043-21	1mm	4.3mmD	\$85.00
6043-22	2mm	4.3mmD	\$85.00
6043-23	3mm	4.3mmD	\$85.00
6050-21	1mm	5.0mmD	\$85.00
6050-22	2mm	5.0mmD	\$85.00
6050-23	3mm	5.0mmD	\$85.00
6060-21	1mm	6.0mmD	\$85.00
6060-22	2mm	6.0mmD	\$85.00
6060-23	3mm	6.0mmD	\$85.00

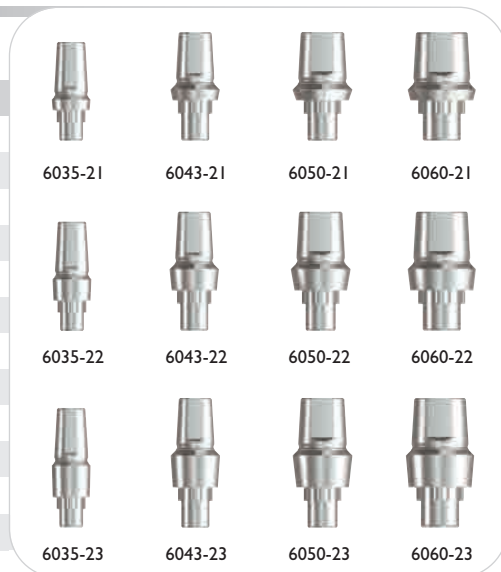

Straight Abutment Transfers

Code	Platform	Price (includes two pieces)
3037-02	3.7mmD	\$30.00
3047-02	4.7mmD	\$30.00
3057-02	5.7mmD	\$30.00
3065-02	6.5mmD	\$30.00



Straight Abutment Comfort Caps

Code	Platform	Price (includes two pieces)
3037-26	3.7mmD	\$30.00
3047-26	4.7mmD	\$30.00
3057-26	5.7mmD	\$30.00
3065-26	6.5mmD	\$30.00



Straight Abutment Analogs

Includes Castable Coping

Code	Platform	Price
3037	3.7mmD	\$35.00
3047	4.7mmD	\$35.00
3057	5.7mmD	\$35.00
3065	6.5mmD	\$35.00



Straight Abutment Castable Copings

Code	Platform	Price (includes two pieces)
3037-28	3.7mmD	\$30.00
3047-28	4.7mmD	\$30.00
3057-28	5.7mmD	\$30.00
3065-28	6.5mmD	\$30.00





Selection Chart: Abutments for Partially Edentulous Restorations

3.5mmD Platform



4.3mmD Platform



5.0mmD Platform



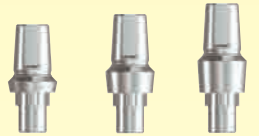
6.0mmD Platform




Straight Abutments - Includes Transfer and Comfort Cap



6035-21 6035-22 6035-23



6043-21 6043-22 6043-23

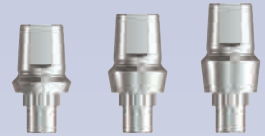

6050-21 6050-22 6050-23

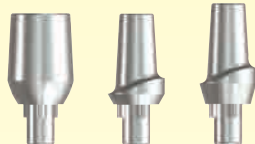


6060-21 6060-22 6060-23


Straight Preparable and Contoured Abutments



6035-30L 6035-31 6035-32



6043-30L 6043-31 6043-32

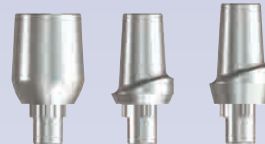


6050-30L 6050-31 6050-32




6060-30L 6060-31 6060-32

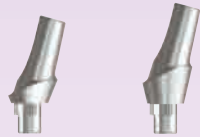
Angled Contoured Abutments




6035-51 6035-52

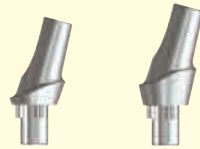

6043-51 6043-52

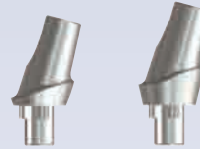

6050-51 6050-52

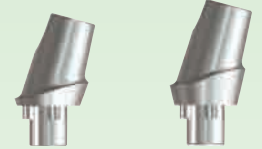

6060-51 6060-52

Straight Contoured Abutments: Tri-Lobe at zero and 30/60/90 Degrees

30/60/90 Degree
6035-32-30
6035-32-60
6035-32-90



30/60/90 Degree
6043-32-30
6043-32-60
6043-32-90



30/60/90 Degree
6050-32-30
6050-32-60
6050-32-90



30/60/90 Degree
6060-32-30
6060-32-60
6060-32-90



Angled Contoured Abutments: Tri-Lobe at zero and 30/60/90 Degrees

30/60/90 Degree
6035-52-30
6035-52-60
6035-52-90



30/60/90 Degree
6043-52-30
6043-52-60
6043-52-90



30/60/90 Degree
6050-52-30
6050-52-60
6050-52-90



Engaging and Non-Engaging Abutments



6035-40 †6035-41 †6035-46 6035-47 †6035-49



6043-40 †6043-41 †6043-46 6043-47 †6043-49



6050-40 †6050-41 †6050-46 6050-47 †6050-49



†6060-46 6060-47 †6060-49

†non-engaging



Healing Collars, Transfers, Analogs & Alignment Tools

Healing Collars

Requires 1.25mmD Hex Tool

Code	Width	Length	Platform	Price
6035-13	4.6mmD	3mmL	3.5mmD	\$30.00
6035-15	4.6mmD	5mmL	3.5mmD	\$30.00
6043-13	5.6mmD	3mmL	4.3mmD	\$30.00
6043-15	5.6mmD	5mmL	4.3mmD	\$30.00
6050-13	6.1mmD	3mmL	5.0mmD	\$30.00
6050-15	6.1mmD	5mmL	5.0mmD	\$30.00
6060-13	6.6mmD	3mmL	6.0mmD	\$30.00
6060-15	6.6mmD	5mmL	6.0mmD	\$30.00



Narrow Healing Collars

Components include Dual-Grip™ feature for use with 1.25mmD Hex Tool or Nobel Unigrip ScrewDriver

Code	Width	Length	Platform	Price
6035-13N	3.5mmD	3mmL	3.5mmD	\$30.00
6035-15N	3.5mmD	5mmL	3.5mmD	\$30.00
6043-13N	4.3mmD	3mmL	4.3mmD	\$30.00
6043-15N	4.3mmD	5mmL	4.3mmD	\$30.00
6050-13N	5.0mmD	3mmL	5.0mmD	\$30.00
6050-15N	5.0mmD	5mmL	5.0mmD	\$30.00
6060-13N	6.0mmD	3mmL	6.0mmD	\$30.00
6060-15N	6.0mmD	5mmL	6.0mmD	\$30.00



Closed-Tray Implant Transfers

Includes Dual-Grip™ Fixation Screw


Code	Platform	Replacement Screw - 2ea	Price
6035-05	3.5mmD	1000-28	\$30.00
6043-05	4.3mmD	1000-27	\$30.00
6050-05	5.0mmD	1000-27	\$30.00
6060-05	6.0mmD	1000-27	\$30.00



Open-Tray Implant Transfers

Includes Long Fixation Screw [Requires 1.25mmD Hex Tool]

Code	Platform	Replacement Screw - 2ea	Price
6035-08	3.5mmD	1000-12	\$30.00
6043-08	4.3mmD	1000-11	\$30.00
6050-08	5.0mmD	1000-11	\$30.00
6060-08	6.0mmD	1000-11	\$30.00

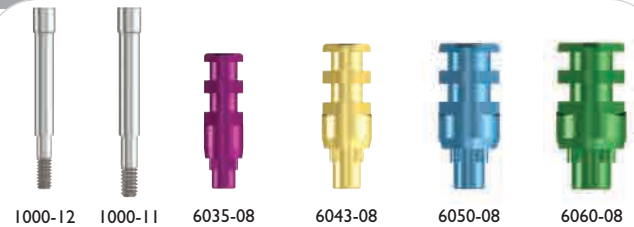


Implant Analogs

Code	Platform	Price
6035-06	3.5mmD	\$20.00
6043-06	4.3mmD	\$20.00
6050-06	5.0mmD	\$20.00
6060-06	6.0mmD	\$20.00



Lab Abutment Alignment Tools

Code	Platform	Price
6035-07	3.5mmD	\$40.00
6043-07	4.3mmD	\$40.00
6050-07	5.0mmD	\$40.00
6060-07	6.0mmD	\$40.00

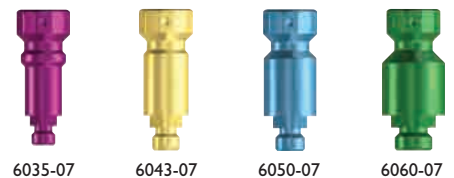



Engaging and Non-Engaging Abutments

Gold Engaging Cast-to Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-40	1mm	3.5mmD	\$100.00
6043-40	1mm	4.3mmD	\$100.00
6050-40	1mm	5.0mmD	\$100.00



6035-40



6043-40



6050-40

Gold Non-Engaging Cast-to Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-41	1mm	3.5mmD	\$100.00
6043-41	1mm	4.3mmD	\$100.00
6050-41	1mm	5.0mmD	\$100.00



6035-41



6043-41



6050-41

Plastic Non-Engaging Castable Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-46	0mm	3.5mmD	\$30.00
6043-46	0mm	4.3mmD	\$30.00
6050-46	0mm	5.0mmD	\$30.00
6060-46	0mm	6.0mmD	\$30.00



6035-46



6043-46



6050-46



6060-46

Plastic Engaging Temporary Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-47	0mm	3.5mmD	\$35.00
6043-47	0mm	4.3mmD	\$35.00
6050-47	0mm	5.0mmD	\$35.00
6060-47	0mm	6.0mmD	\$35.00



6035-47



6043-47



6050-47



6060-47

Titanium Non-Engaging Temporary Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-49	0mm	3.5mmD	\$35.00
6043-49	0mm	4.3mmD	\$35.00
6050-49	0mm	5.0mmD	\$35.00
6060-49	0mm	6.0mmD	\$35.00



6035-49



6043-49



6050-49



6060-49

Straight Wide Lab Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-30L	0mm	3.5mmD	\$75.00
6043-30L	0mm	4.3mmD	\$75.00
6050-30L	0mm	5.0mmD	\$75.00
6060-30L	0mm	6.0mmD	\$75.00



6035-30L



6043-30L



6050-30L



6060-30L

Straight Contoured Abutments

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-31	1mm	3.5mmD	\$85.00
6035-32	2mm	3.5mmD	\$85.00
6043-31	1mm	4.3mmD	\$85.00
6043-32	2mm	4.3mmD	\$85.00
6050-31	1mm	5.0mmD	\$85.00
6050-32	2mm	5.0mmD	\$85.00
6060-31	1mm	6.0mmD	\$85.00
6060-32	2mm	6.0mmD	\$85.00



6035-31



6043-31



6050-31



6060-31



6035-32



6043-32



6050-32



6060-32




Engaging and Non-Engaging Abutments - Continued

Angled Contoured Abutments

Includes Dual-Grip™ Fixation Screw


Code	Collar Heights	Platform	Price
6035-51	.75mm	3.5mmD	\$85.00
6035-52	2mm	3.5mmD	\$85.00
6043-51	.75mm	4.3mmD	\$85.00
6043-52	2mm	4.3mmD	\$85.00
6050-51	.75mm	5.0mmD	\$85.00
6050-52	2mm	5.0mmD	\$85.00
6060-51	.75mm	6.0mmD	\$85.00
6060-52	2mm	6.0mmD	\$85.00

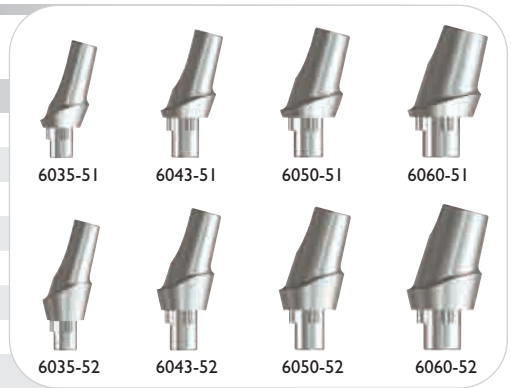

Angled Contoured Abutments with 30, 60 & 90 Degree offset Tri-Lobe

Includes Dual-Grip™ Fixation Screw


Code	Collar Heights	Platform	Price
6035-52-30	2mm	3.5mmD	\$85.00
6035-52-60	2mm	3.5mmD	\$85.00
6035-52-90	2mm	3.5mmD	\$85.00
6043-52-30	2mm	4.3mmD	\$85.00
6043-52-60	2mm	4.3mmD	\$85.00
6043-52-90	2mm	4.3mmD	\$85.00
6050-52-30	2mm	5.0mmD	\$85.00
6050-52-60	2mm	5.0mmD	\$85.00
6050-52-90	2mm	5.0mmD	\$85.00

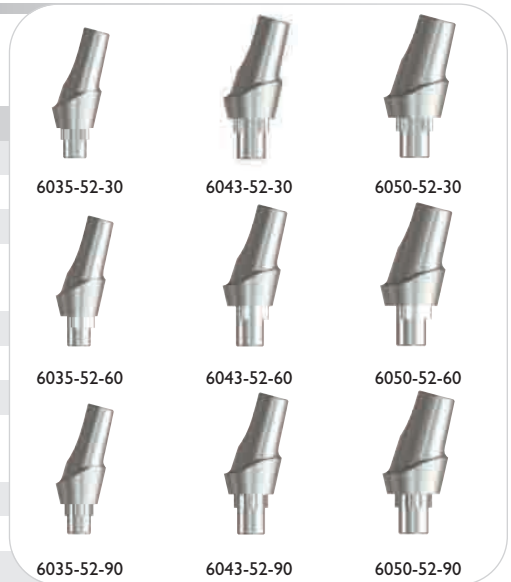

Straight Contoured Abutments with 30, 60 & 90 Degree offset Tri-Lobe

Includes Dual-Grip™ Fixation Screw

Code	Collar Heights	Platform	Price
6035-32-30	2mm	3.5mmD	\$85.00
6035-32-60	2mm	3.5mmD	\$85.00
6035-32-90	2mm	3.5mmD	\$85.00
6043-32-30	2mm	4.3mmD	\$85.00
6043-32-60	2mm	4.3mmD	\$85.00
6043-32-90	2mm	4.3mmD	\$85.00
6050-32-30	2mm	5.0mmD	\$85.00
6050-32-60	2mm	5.0mmD	\$85.00
6050-32-90	2mm	5.0mmD	\$85.00
6060-32-30	2mm	6.0mmD	\$85.00
6060-32-60	2mm	6.0mmD	\$85.00
6060-32-90	2mm	6.0mmD	\$85.00

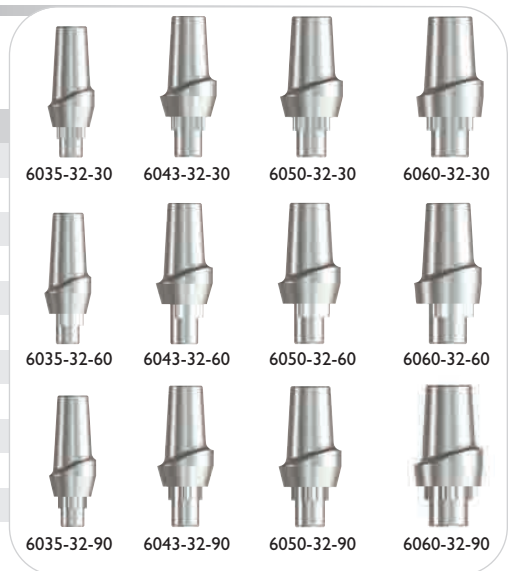

Replacement Fixation Screws

Components include Dual-Grip™ feature for use with 1.25mmD Hex Tool or Nobel Unigrip ScrewDriver



Code	Platform	Price (includes two pieces)
1000-23	3.5mmD only (short)	\$30.00
1000-28	3.5mmD only (long)	\$30.00
1000-22	4.3, 5.0 and 6.0mmD (short)	\$30.00
1000-27	4.3, 5.0 and 6.0mmD (long)	\$30.00





Orienting the 15° Angled Abutment

The flat side of the Fixture Mount/Transfer [Fig. 1] and the angle of the 15° abutment [Fig. 2] are manufactured to align with the flat of the implant's internal tri-lobe.



To establish the optimum angle of the abutment, index the flat surface of the Fixture Mount/Transfer either toward or away from the desired direction of abutment angle at the time of implant placement.








Fig. 1
Fixture Mount/Transfer


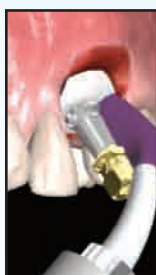












Fig. 2 The two-piece, 15° Angled Abutment consists of an angled body and a fixation screw

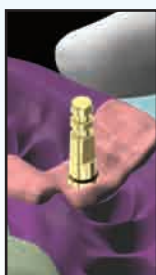

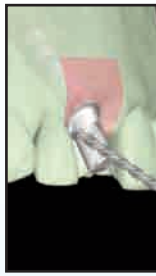


Option 1: Attaching an Angled Abutment and **chairside procedures** for preparation

	Remove the Healing Collar with a 1.25mmD Hex Tool		Attach the abutment with the 1.25mmD Hex Tool		Mark the required modifications to achieve appropriate clearance		Detach the abutment with the 1.25mmD Hex Tool		Attach the abutment to the corresponding Implant Analog. Reduce the abutment to achieve appropriate clearance
--	---	---	---	---	--	--	---	---	---

Option 2: Stage-One transfer followed by **laboratory procedures** for abutment preparation and provisionalization


	During implant placement, index the flat of the Fixture-Mount		Block out the top of the Fixture-Mount and make a conventional transfer impression		Attach the Cover Screw and submerge implant		Attach the Fixture-Mount to the Implant Analog and seat into the negative of the impression		Pour the working cast in dental stone with soft tissue replication material		Remove the Fixture-Mount from the working cast with the 1.25mmD Hex Tool
	Attach the Angled Abutment to the Implant Analog in the working cast with a 1.25mmD Hex Tool		Modify angle of abutment and refine margins as required		Fabricate the temporary crown		At Stage-2 Surgery, attach the prepared abutment, tighten to 30 Ncm with a Torque Wrench		Occlude screw access hole and confirm contour of restoration		Cement the provisional restoration

Common procedures for abutment placement, margin finessing, impressing and placement of restoration


	Tighten to 30 Ncm with a Torque Wrench		Finesse the margins with a high-speed bur and copious irrigation prior to final impressing		Make an impression following standard clinical procedures		Attach provisional restoration and cement onto the prepared abutment		Follow standard laboratory procedures, fabricate the final restoration and cement onto the prepared abutment
--	--	---	--	---	---	--	--	---	--



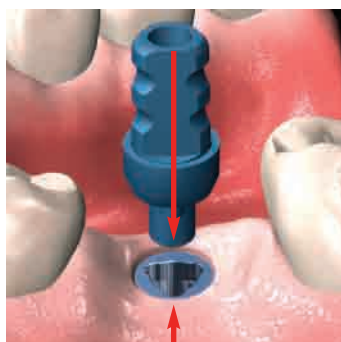Contoured Abutments with 0, 30, 60 or 90 Degree Tri-Lobe Offset



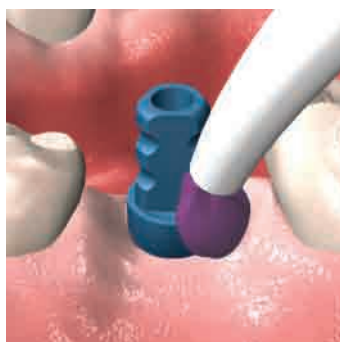
Implant with tri-lobe properly aligned to labial. Rotate transfer and seat mating tri-lobes.




Make transfer impression using vinyl polysiloxane material.



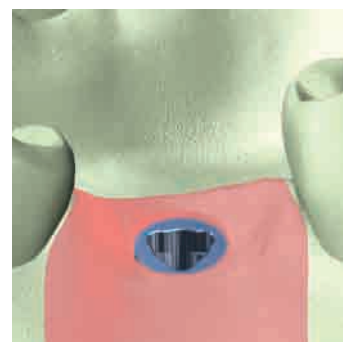
Attach transfer to matching colored implant analog. Insert transfer and pour impression.



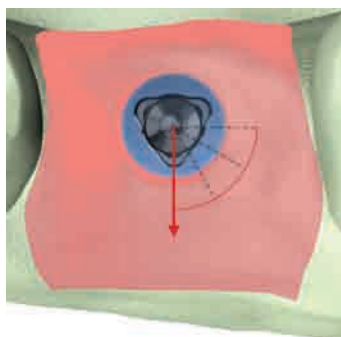
Evaluate rotational position of tri-lobe for selection of straight and angled contoured abutments



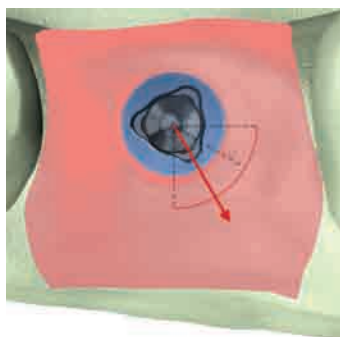
Tri-Lobe aligned perpendicular to labial. (0 degrees offset)



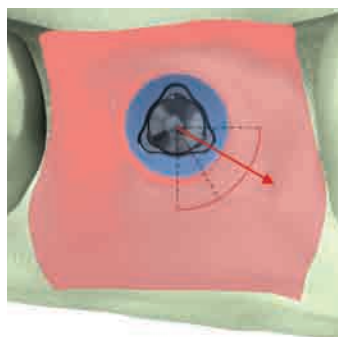
Tri-Lobe aligned 30 degrees counter-clockwise to line perpendicular to labial.



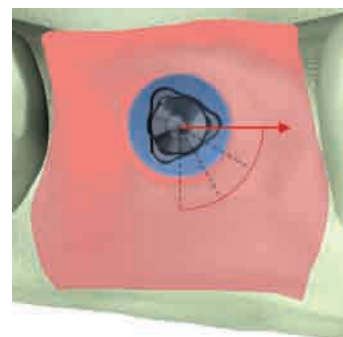
Tri-Lobe aligned 60 degrees counter-clockwise to line perpendicular to labial.



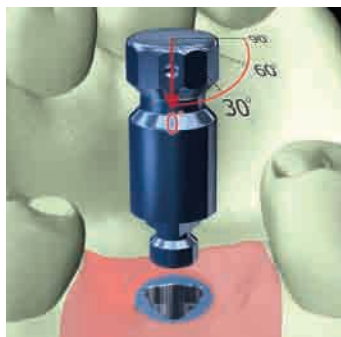
Tri-Lobe aligned 90 degrees counter-clockwise to line perpendicular to labial.



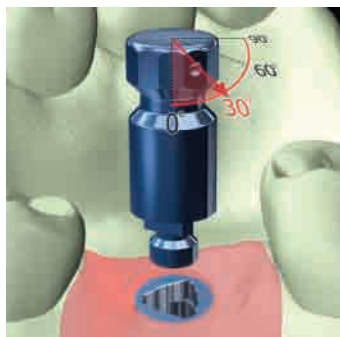
Dimple of Alignment Tool is perpendicular to labial: 0 degree offset of implant tri-lobe.



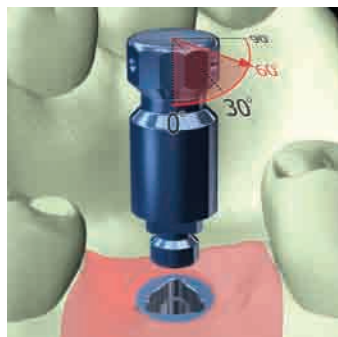
Dimple on the Tool is 1 facet (30 degrees) counter-clockwise to line perpendicular to labial.



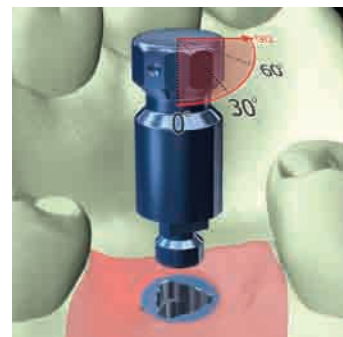
Dimple on the Tool is 2 facets (60 degrees) counter-clockwise to line perpendicular to labial.



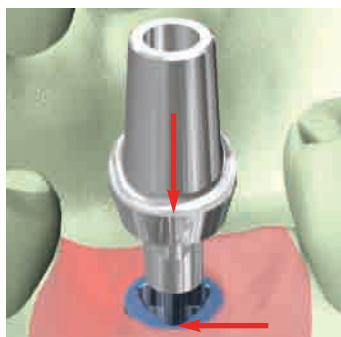
Dimple on the tool is 3 facets (90 degrees) counter-clockwise to line perpendicular to labial.



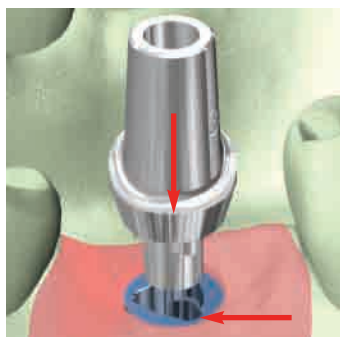
Use 0 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



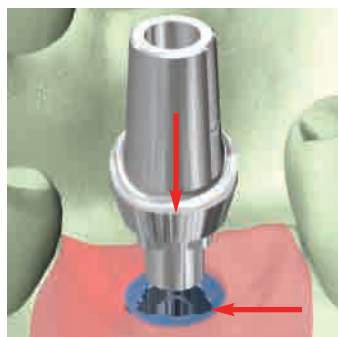Use 30 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



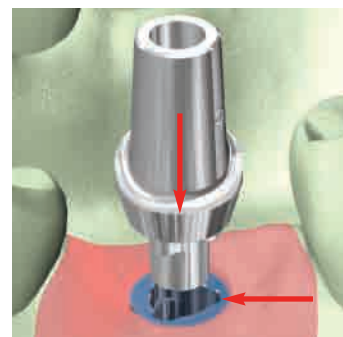Use 60 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



Use 90 degree tri-lobe offset Contoured Abutment when tool is inserted as above.



Selection Chart: Overdenture & Screw-Retained Restorations

	3.5mmD Platform	4.3mmD Platform	5.0mmD Platform	6.0mmD Platform				
	Overdenture Abutments (*Includes Transfer and Comfort Cap)							
Overdenture Abutments	 6035-61 6035-62 6035-64	 6035-64-30 6035-64-15	 6043-61 6043-62 6043-64	 6043-64-30 6043-64-15	 6050-61 6050-62 6050-64	 6050-64-30 6050-64-15	 6060-61 6060-62 6060-64	
	Overdenture Abutment Transfers, Analog & Comfort Cap							
OD Abutment Comfort Cap Transfers and Analog	 *4000-13	 *1000-68	 1000-61	 1000-62	 1000-63	 1000-66	 1000-67	 1000-69
	Overdenture Copings							
Overdenture Abutment Components	 1000-64 Gold	 1000-65 Titanium	 1000-66 Plastic	 1000-67 Gold/Plastic	 1000-69 Titanium			
	Ball Abutments (**Includes Transfer and Comfort Cap)							
Ball Abutments	 6035-72 6035-74	 6043-72 6043-74	 6050-72 6050-74	 6060-72 6060-74				
	Ball Attachments							
Ball Abutment Components	 **1000-87	 **1000-84	 1000-85	 1000-81	 1000-82	 1000-88	 1000-83	 1000-86
	Locator® Abutments							
Locator® Abutments	 6035-90 6035-91 6035-92 6035-93 6035-94	 6043-90 6043-91 6043-92 6043-93 6043-94	 6050-91 6050-92 6050-93 6050-94					
	Locator® Attachments							
Locator® Abutment Components	 1000-94	 1000-95	 1000-92B 1000-92C	 1000-92G 1000-92P	 1000-92R	 1000-98	 1000-93	 1000-96




Abutments for Ball, Bar & Locator® Overdentures

Overdenture Abutments

Includes Transfer & Comfort Cap
[Requires 1.25mmD Hex Tool]

Code	Collar Heights	Platform	Price
6035-61	1mm	3.5mmD	\$85.00
6035-62	2mm	3.5mmD	\$85.00
6035-64	4mm	3.5mmD	\$85.00
6043-61	1mm	4.3mmD	\$85.00
6043-62	2mm	4.3mmD	\$85.00
6043-64	4mm	4.3mmD	\$85.00
6050-61	1mm	5.0mmD	\$85.00
6050-62	2mm	5.0mmD	\$85.00
6050-64	4mm	5.0mmD	\$85.00
6060-61	1mm	6.0mmD	\$85.00
6060-62	2mm	6.0mmD	\$85.00
6060-64	4mm	6.0mmD	\$85.00

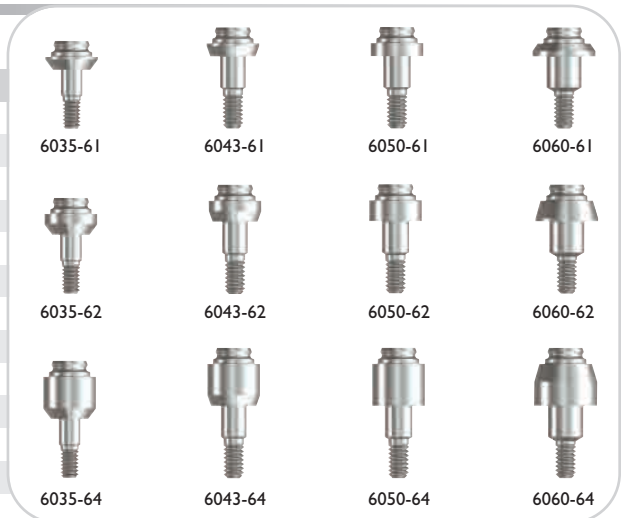

Angled Overdenture Abutments

Includes Dual-Grip™ Abutment Screw, Transfer & Comfort Cap [Screw Top requires 2.5mmD Hex Tool]


Code	Description	Platform	Price
6035-64-15	15° Angled Abutment	3.5mmD	\$100.00
6043-64-15	15° Angled Abutment	4.3mmD	\$100.00
6050-64-15	15° Angled Abutment	5.0mmD	\$100.00
6035-64-30	30° Angled Abutment	3.5mmD	\$100.00
6043-64-30	30° Angled Abutment	4.3mmD	\$100.00
6050-64-30	30° Angled Abutment	5.0mmD	\$100.00
1000-60	Replacement Screw Top for Angled Base		\$35.00



Ball Abutments


Includes Transfer & Comfort Cap
[Requires 1.25mmD Hex Tool]

Code	Collar Heights	Platform	Price
6035-72	2mm	3.5mmD	\$85.00
6035-74	4mm	3.5mmD	\$85.00
6043-72	2mm	4.3mmD	\$85.00
6043-74	4mm	4.3mmD	\$85.00
6050-72	2mm	5.0mmD	\$85.00
6050-74	4mm	5.0mmD	\$85.00
6060-72	2mm	6.0mmD	\$85.00
6060-74	4mm	6.0mmD	\$85.00
1000-70	Replacement Ball Top for Angled Base		\$35.00

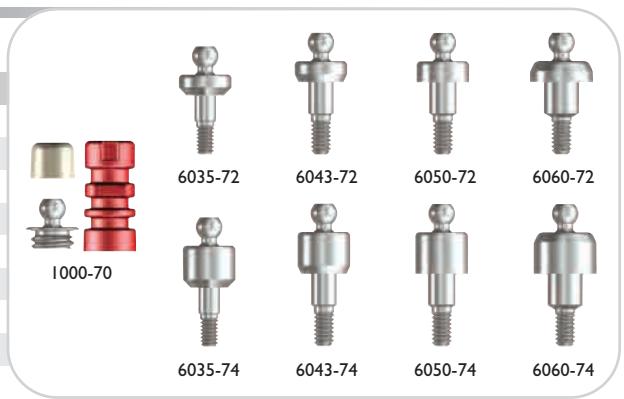

RePlant Locator® Abutments

Code	Collar Heights	Platform	Price
6035-90	0mm	3.5mmD	\$110.00
6035-91	1mm	3.5mmD	\$110.00
6035-92	2mm	3.5mmD	\$110.00
6035-93	3mm	3.5mmD	\$110.00
6035-94	4mm	3.5mmD	\$110.00
6043-90	0mm	4.3mmD	\$110.00
6043-91	1mm	4.3mmD	\$110.00
6043-92	2mm	4.3mmD	\$110.00
6043-93	3mm	4.3mmD	\$110.00
6043-94	4mm	4.3mmD	\$110.00
6050-91	1mm	5.0mmD	\$110.00
6050-92	2mm	5.0mmD	\$110.00
6050-93	3mm	5.0mmD	\$110.00
6050-94	4mm	5.0mmD	\$110.00

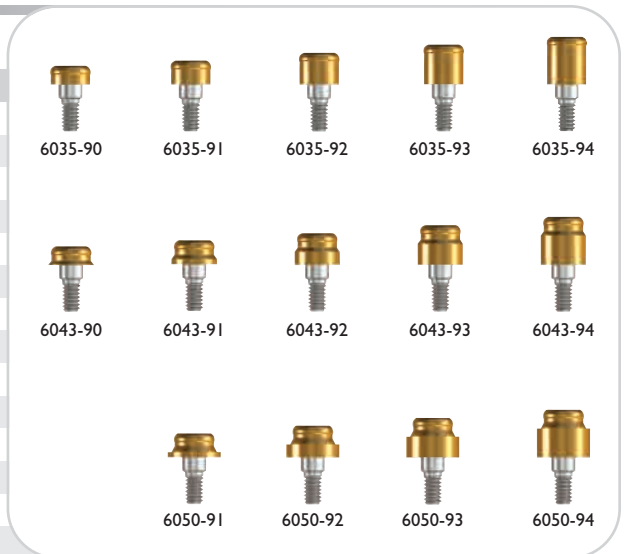




Ancillary Components for Overdenture & Screw-Retained Restorations

Overdenture Abutment Comfort Cap, Transfers & Analog Screws require 1.25mmD Hex Tool

Code	Platform	Description	Price
1000-61	5.0mmD	Aluminum Analog	\$20.00
1000-62	5.0mmD	Open-Tray Transfer plus Fixation Screw	\$30.00
1000-63	5.0mmD	Closed-Tray Transfer plus Fixation Screw	\$30.00
1000-68	5.0mmD	Plastic Abutment Transfer - 2ea	\$30.00
4000-13	5.0mmD	Plastic Comfort Cap - 2ea	\$30.00
1000-14	-	Replacement Open-Tray Fixation Screw - 2ea	\$30.00
1000-15	-	Replacement Closed-Tray Fixation Screw - 2ea	\$30.00

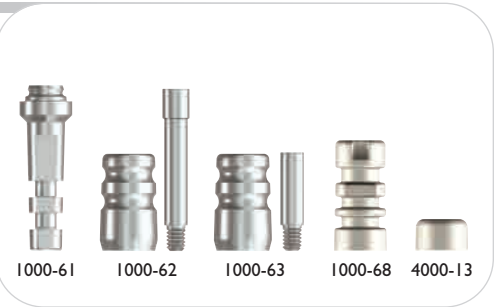


Overdenture Abutment Components Includes Fixation Screw [1000-05], optional Dual-Grip™ Fixation Screw [1000-25]

Code	Platform	Description	Price
1000-64	5.0mmD	Gold Coping	\$80.00
1000-65	5.0mmD	Titanium Coping	\$40.00
1000-66	5.0mmD	Plastic Coping	\$35.00
1000-67	5.0mmD	Gold & Plastic Coping	\$80.00
1000-69	5.0mmD	Titanium Temporary Coping	\$35.00
1000-05	-	Coping Fixation Screw - 2ea	\$30.00
1000-25	-	Dual-Grip™ Coping Fixation Screw - 2ea	\$30.00



Ball Abutment Components

Code	Description	Price
1000-81	Ball Abutment Housing - 2ea	\$30.00
1000-82	Nylon Liner - 2ea	\$15.00
1000-83	Nylon Liner Transfer - 2ea	\$15.00
1000-84	Ball Abutment Transfer - 2ea	\$30.00
1000-85	Ball Abutment Analog	\$20.00
1000-86	Ball Abutment Tools	\$65.00
1000-87	Ball Abutment Comfort Cap - 2ea	\$30.00
1000-88	Ball Abutment Housing & Nylon Liner	\$22.00

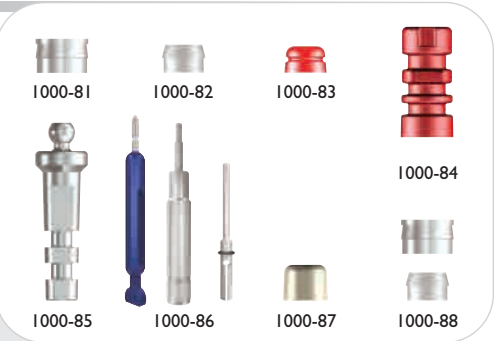

Locator® Abutment Transfer, Analogs and Alignment Tools

Code	Description	Price
1000-94	Locator® Impression Coping - 2ea	\$32.00
1000-95	Locator® Abutment Analog 4mmD (for 3.5 & 4.3mmD) - 2ea	\$32.00
1000-95W	Locator® Abutment Analog 5mmD (for 5.0mmD) - 2ea	\$32.00
1000-97	Locator® Abutment Parallel Post - 4ea	\$16.00
1000-11	Guide Screw for 4.3, 5.0 and 6.0mmD implant platform - 2ea	\$30.00
1000-12	Guide Screw for 3.5mmD implant platform only - 2ea	\$30.00
1000-99	Locator® Abutment Angle Measure Guide	\$21.00



Locator® Abutment Components

Code	Retention Force	Description	Price
1000-92C	5.0lbs	Replacement Male Clear - 4ea	\$21.00
1000-92P	3.0lbs	Light Replacement Male Pink - 4ea	\$21.00
1000-92B	1.5lbs	Extra Light Replacement Male Blue - 4ea	\$21.00
1000-92G	4.0lbs	Extended Range Male Green - 4ea	\$21.00
1000-92R	1.0lb	Extra Light Extended Range Male Red - 4ea	\$21.00
1000-93	-	Black Processing Male - 4ea	\$21.00
1000-98	Assorted	Abutment Processing Cap/Male Kit - 2ea	\$58.00



Locator® Instrumentation

Code	Description	Price
LDT	Locator® Abutment Driver Tool - Long	\$44.00
LDTs	Locator® Abutment Driver Tool - Short	\$44.00
LDD	Locator® Abutment Driver Drill - Long	\$44.00
LDdS	Locator® Abutment Driver Drill - Short	\$44.00
1000-96	Core Tool - Male Removal & Seating Tool and Abutment Driver	\$75.00





Tapered Hex Tools

Code	Description	Price
THT1.25	1.25mmD Tapered Hex Tool Long	\$40.00
THT1.25S	1.25mmD Tapered Hex Tool Short	\$40.00



THT1.25



THT1.25S

Straight Tools

Square connection for use with Ratchet or Torque Wrench

Code	Description	Price
HT1.25	1.25mmD Straight Hex Tool Long	\$40.00
HT1.25S	1.25mmD Straight Hex Tool Short	\$40.00
HT2.5	2.5mmD Hex Tool Long (with retentive O-ring)	\$40.00
HT2.5S	2.5mmD Hex Tool Short (with retentive O-ring)	\$40.00



HT1.25



HT2.5



HT1.25S



HT2.5S

Hex Drills

Latch head connection for use with Surgical Motor

Code	Description	Price
HD1.25	1.25mmD Hex Drill	\$40.00
HD2.5	2.5mmD Hex Drill (with retentive O-ring)	\$40.00



HD1.25



HD2.5

Ancillary Instruments

Code	Description	Price
TW30	30Ncm Torque Wrench	\$125.00
TW35	35Ncm Torque Wrench	\$125.00
DE	Drill Extender	\$40.00
TWM	Multi-Setting Torque Wrench (15, 20, 25, 30, 35, 40, 45, 50 & 60Ncm)	\$395.00
ABH	Analog Holder	\$85.00



TW30



TW35



DE



TWM



ABH



simply smarter.
 Implant Direct International Inc.
 8840 Russell Road,
 Las Vegas, NV 89148
 Manufacturing/Sales: 1-818-444-3333
 www.implantdirect.com



"Order Online for High Quality Products at factory-direct prices"

Gerald Niznick, DMD, MSD, Founder and President of Implant Direct

ReActive™ Tri-Lobe Implant
 Platform Shifting, Micro-threads,
 Deeper Threads, Round Apex
 All-in-one Packaging: Cover Screw,
 Transfer & Straight Abutment



Surgically and/or Prosthetically Compatible Implants
 Select Type of Threads and Abutment

RePlant with Temporary Titanium Abutment

RePlus with Straight "Snappy" Abutment

ReActive with Straight Preparable Abutment

RePlant®*
 Body Dimensions:
 Matches Replace†

RePlus®
 Body Taper:
 For Bone Expansion

ReActive™
 Body Threads:
 For Increased Surface



*RePlant implant surgically compatible with Nobel Replace® Drills

†Nobel Replace® is a registered Trademark of Nobel Biocare

**Nobel Active® is a registered Trademark of Nobel Biocare

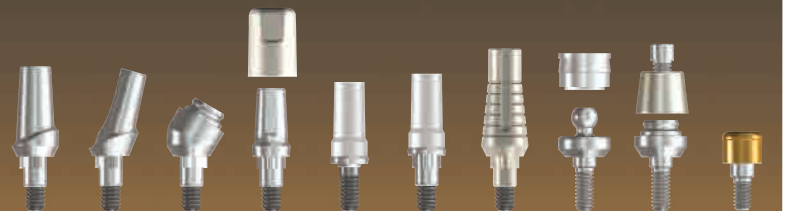
Free Implants for Nobel Biocare Customers
 3 Surgically Compatible, RePlant® Implants with Titanium Transfer/Temporary Abutments plus 2 Contoured Abutments

US List Price for 3 Implants, Abutments and Prosthetic Components:
RePlant = \$620
NobelBioCare = \$1907

8mm/3.5mm RePlant 10mm/3.5mm RePlant 13mm/3.5mm RePlant 2 Contoured Abutments with Fixation Screws

FREE Introductory Offers at
www.implantdirect.com

Terms and Conditions Apply



Compare US List Price	Straight Contoured	15° Angled Contoured	30° Screw Receiving	Straight Snap-On	Gold/ w/o Tri-Lobe	Gold/ with Tri-Lobe	Temporary Abutment Inc. Screw	Ball Attachment	Screw Receiving	Locator® Abutment
Implant Direct	\$85	\$85	\$100	\$85	\$100	\$100	\$35	\$107	\$85	\$110
Nobel Biocare	\$170	\$225	\$269	\$201	\$216	\$216	\$72	\$217	\$229	\$110
Cost Savings	50%	62%	63%	58%	54%	54%	51%	51%	63%	0%

List prices as of 2/23/09

* Stig Hansson. The Implant Neck: Smooth or Provided with Retention Elements. A Biomechanical Approach. Clin Oral Impl Res 1999;10:394-405

NobelBiocare™ and NobelReplace™ are registered Trademarks of the Nobel BioCare Group