



# Zirkonzahn®

*Human Zirconium Technology*

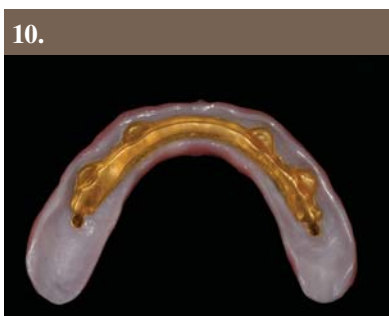

**ALL ZIRCONIA**

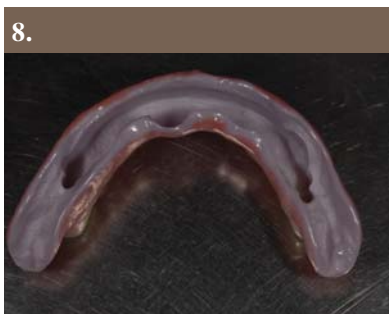

1. Upper and lower primary milled zirconia bars

2. Secondary part produced with galvano technique

3. Mock-up resin frame work

4. Copy-mill resin frames in zirconia

5. Zirconia frame coloured in tissue- and dentin shades

6. Tissue-coloured veneer porcelain and thin layer of dentin-coloured ceramic applied

7. Secondary galvano parts in situ

8. Milled and veneered zirconia frame ready for cementation

9. Connect parts with anaerobic cement

10. Galvano part and bridge permanently joined

11.



12.



11. Bridge glides smoothly over milled bar

12. Bridge seated

13.



14.



13. Case complete

14. Case complete

15.



16.



15. Screw-retained bars


16. Screw-retained bars



Issued case in situ

Zilio Aldo, Venice

# Zirkonzahn®

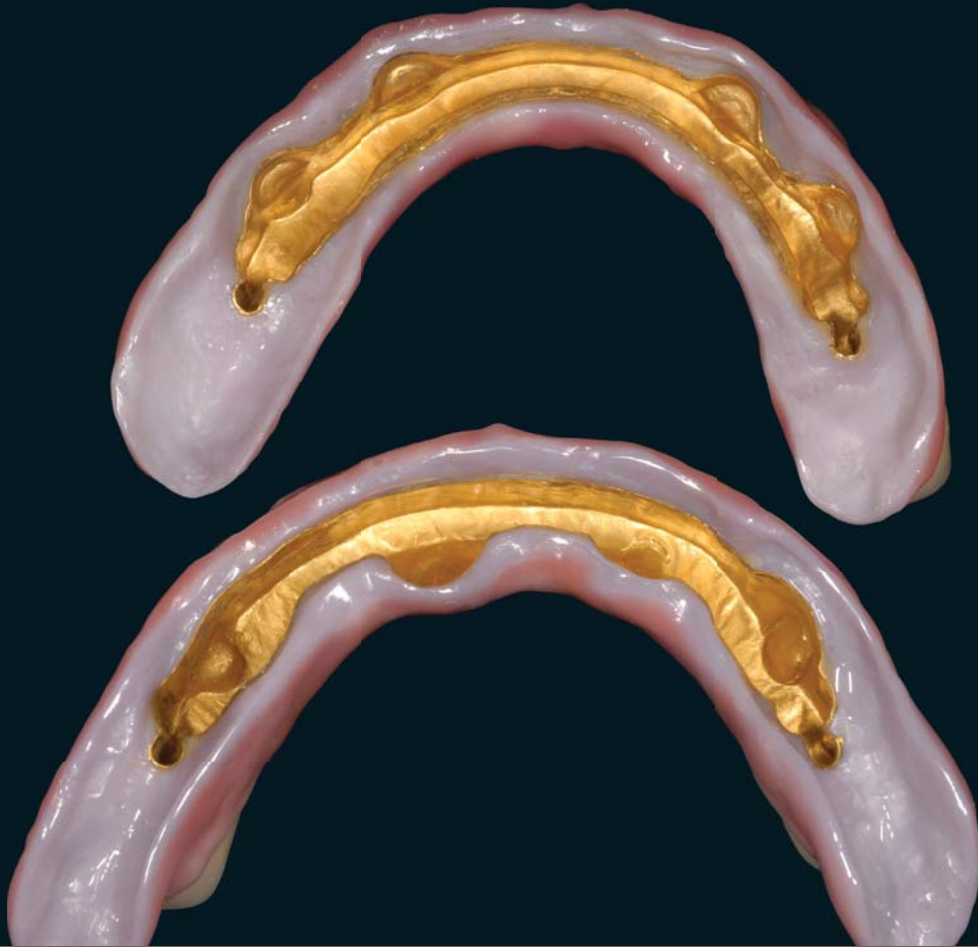

## HUMAN ZIRCONIUM TECHNOLOGY

ZIRKONZAHN USA Inc. - Tel 1 800 989 8931 - Tel 678 441 9419 - [infousa@zirkonzahn.com](mailto:infousa@zirkonzahn.com)

ZIRKONZAHN World Wide - Tel +39 0474 066 680 - [info@zirkonzahn.com](mailto:info@zirkonzahn.com) - [www.zirkonzahn.com](http://www.zirkonzahn.com)



*milled with Zirkonzahn Milling System*



WEAA5401=