


# Professional abrasives meeting the highest demands

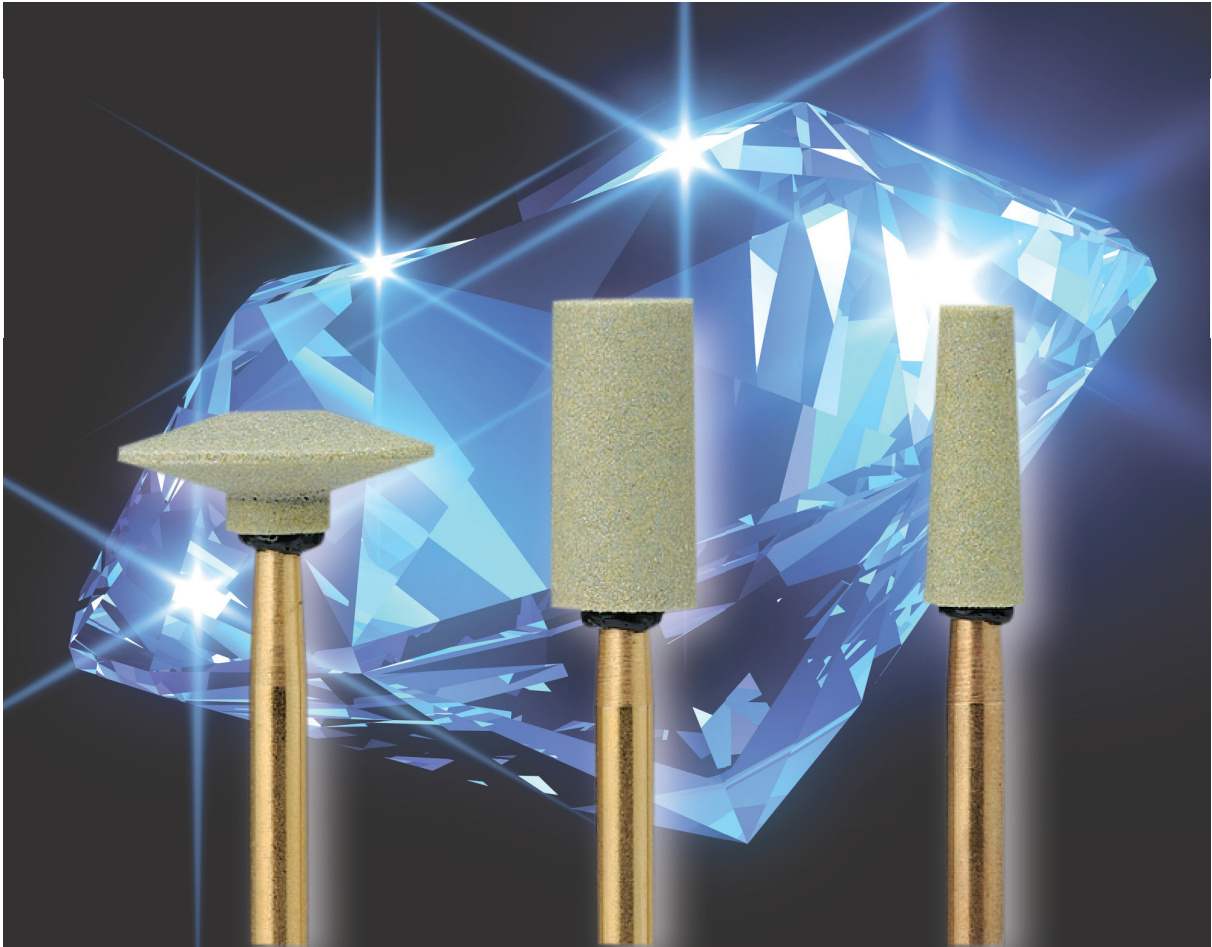

## Dura-Green® DIA

The contouring and finishing of all-ceramic and zirconia-based restorations has now become much easier, quicker and safer.

Using diamond particles in a special binder, **Dura-Green DIA** provide outstanding abrasive properties and high efficiency in a wide range of applications.

**Dura-Green DIA** abrasives will allow you to adjust all types of ceramic surfaces, fine anatomical details or fragile, sensitive marginal areas at a low contact pressure – minimising the risk of chipping or micro-cracks.



Proven products for better dentistry

# Dura-Green® DIA

EFFICIENT, VERSATILE, GENTLE ON THE MATERIAL



**CY5: Surface adjustment**  
Excellent cutting performance offers convenient surface adjustment with efficient material reduction.



**KN7: Framework adjustment**  
Gentle on the material when making interdental or marginal corrections and wall thickness adjustments.



**TC4: Detailed finishing**  
Whether dressed or in its original shape, TC4 helps to efficiently create surface details and make occlusal adjustments.




For cleaning and shape adjustment, press Dura-Green DIA against the Shofu Adjuster (aluminium oxide) using a grinding motion.

## Advantages at a glance:

- Efficient and gentle on the material, thanks to Dura-Green DIA abrasives being impregnated with diamonds throughout
- Minimised risk of micro-cracking
- Low contact pressure
- Low vibration
- Minimal heat generation
- High edge stability

The diamond-impregnated microstructure of Dura-Green DIA abrasives ensures a more homogeneous surface pattern and a longer service life, as compared to silicon carbide instruments:

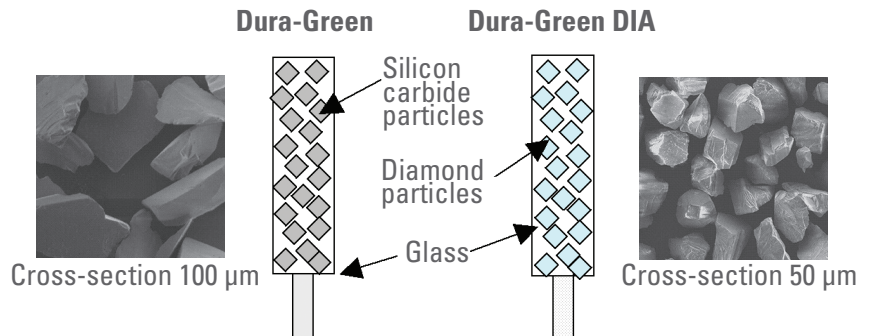

**Ideal for use with:** all-zirconia or zirconia-based restorations, lithium disilicate, aluminium oxide, pressed & layered ceramics and metal surfaces

## SHAPES

Maximum speed:  
30,000 rpm

Shape	CY5	FL3	IC5	KN7	RD3	TC2	TC4	WH6
ISO No.	056	031	070	125	035	025	040	125
HP	0156	0158	0162	0154	0159	0161	0155	0157
CA		0163			0164	0165		

Box of 2 pcs.: CY5, KN7, WH6 / box of 3 pcs.: FL3, IC5, RD3, TC2, TC4

PN	Item	Contents
0153	Kit HP: CY5, KN7, TC4	1 pce. each
0166	Kit CA: FL3, RD3, TC2	1 pce. each
0504	Shofu Adjuster	1 pce.



**SHOFU DENTAL GMBH**  
Am Brüll 17 · D-40878 Ratingen  
Telefon 0 21 02/86 64-0 · Fax 0 21 02/86 64-65  
E-mail: info@shofu.de · www.shofu.de