


Tribos 501

Tribos 501 su visokokvalitetni zubi razvijeni u dijagonalnoj mreži PMMA, proizvedeni u novom kondenzacionom postupku „Mega-Press-Inject“ za optimizaciju gustine i tvrdoće površine zuba.

Ova proizvodnja u 3 sloja pruža isti kvalitet u svim područjima zuba.

Ovo posebno utiče na trajnost i osetljivost proteza na plak. Korisnik će osetiti preciznost oblika i tvrdoće Tribosa 501 kada ga prvi put koristi, što je uporedivo s preciznošću CAD / CAM tehnologije.

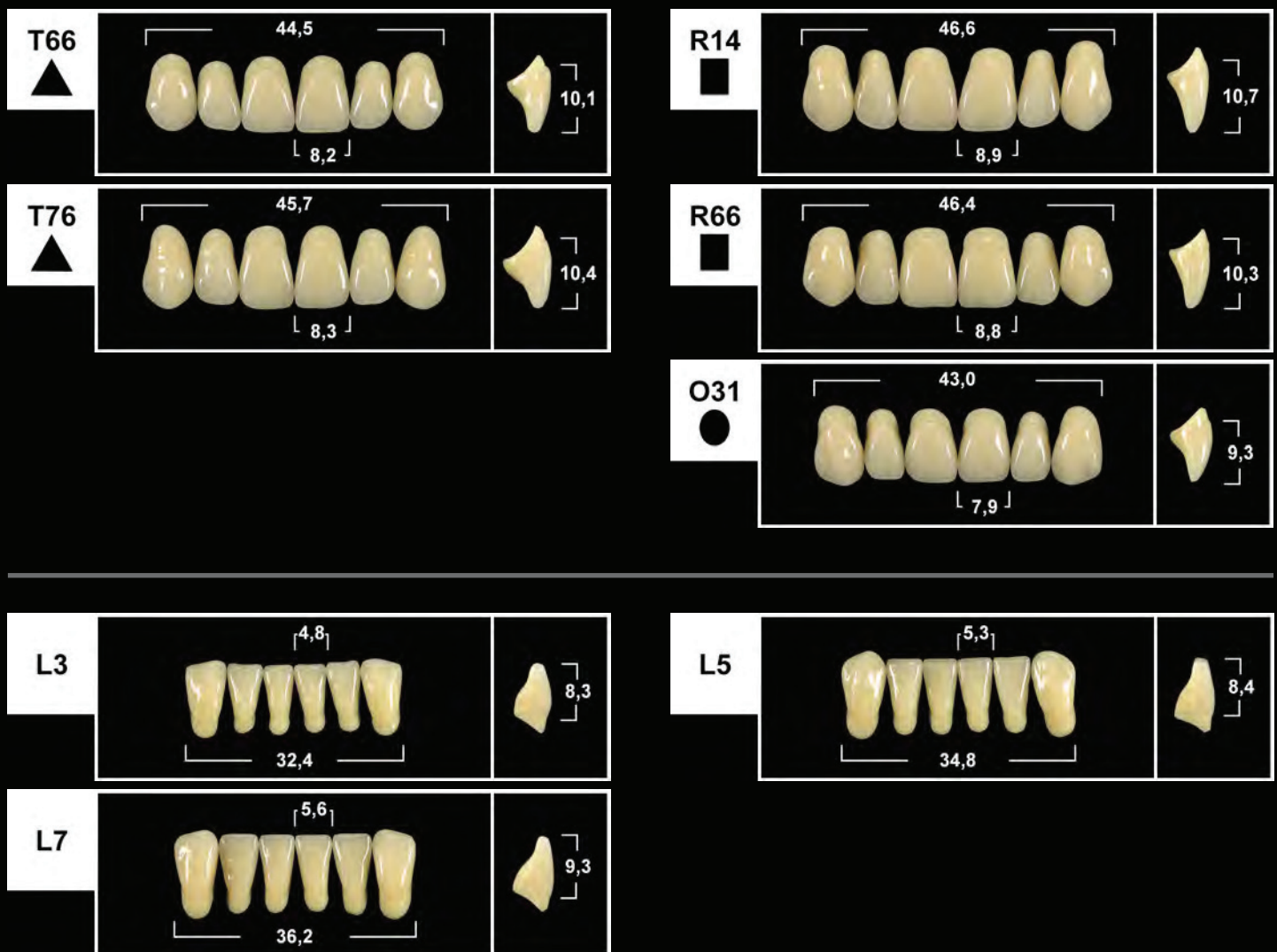

Da bismo koncept NŠP predstavili što bolje, pripremili smo za Vas slučajeve zuba i kartice oblika, pored uputstava o tome kako da koristite ceo koncept. Pojedinačne komplete možete pogledati i zasebno na: vww.gebdi-dental.com i poručiti preko MIKODENTAL Dental Depoa.

Za tehničku pomoć, molimo Vas da kontaktirate našeg partnera: MIKODENTAL, gde će Vam sve biti objašnjeno od strane stručnih lica.

















TRIBOS 501- PREMIJUM

Visokokvalitetni troslojni zubi sa preciznošću CAD/CAM tehnologije!
Morfologija rezultat 3-godišnje studije.

CELA GARNITURA 50€



15€
-40% 25€

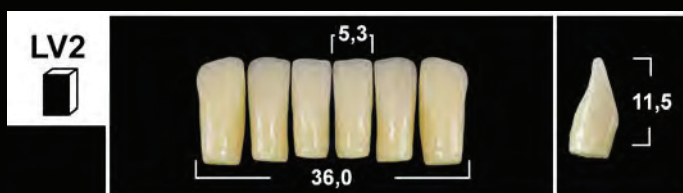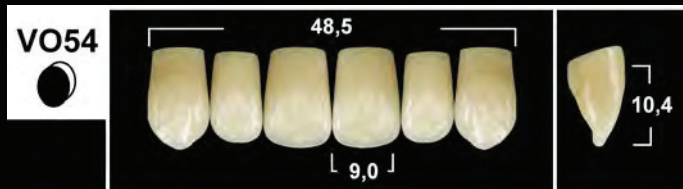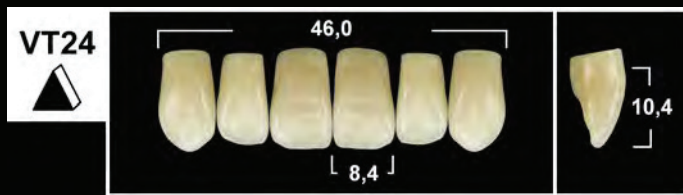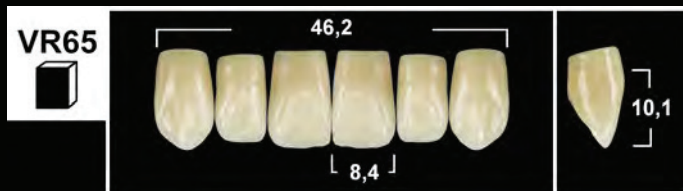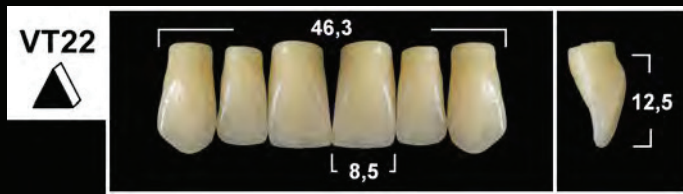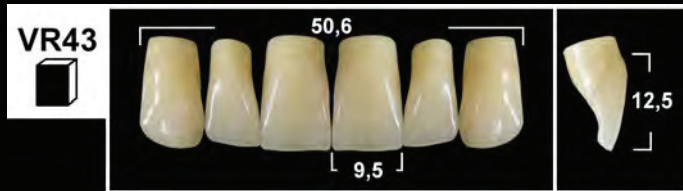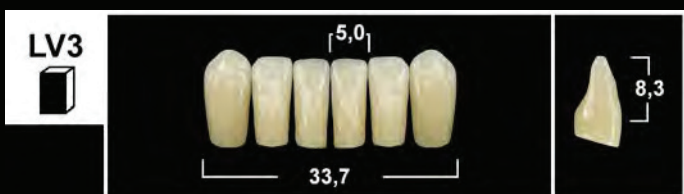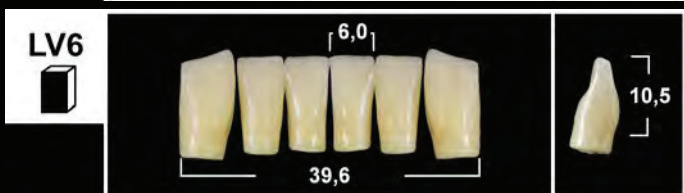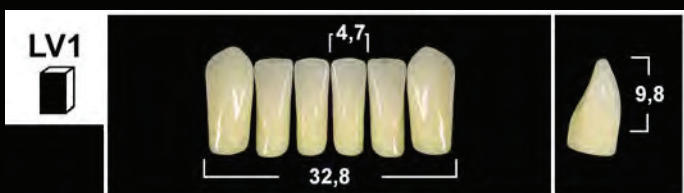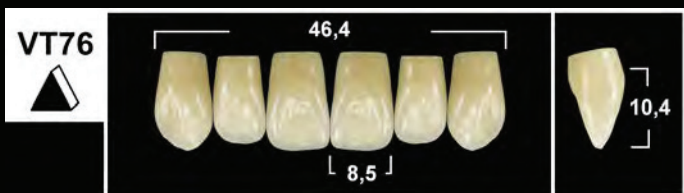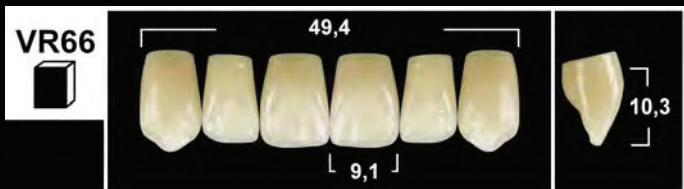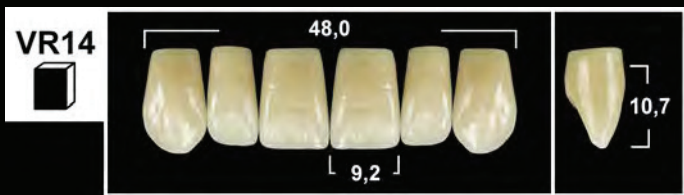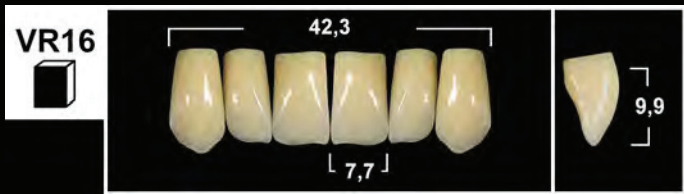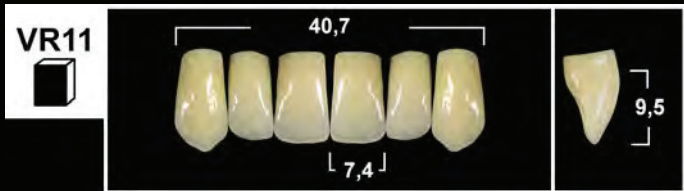

Primarna morfologija - Sekundarna morfologija - Sekundarna morfologija abrazivna

Primarna morfologija



Primarna morfologija (PM / PS) odgovara morfologiji mladog zuba koji nema abraziju, niti kontaktne površine. Osnova zuba je smanjena što omogućava da se koristi u kompletnoj protezi, kao i implantatima.

Montažni sto 140
Sagitalni nagib, bočni Condilar (SKN) FH 45° - CE 30°
Bennetni ugao 10° ISS 0,5 mm
Laterotrusion LRT -20°
Retrusion RT 0,5 mm



Sekundarna morfologija



Sekundarna morfologija (SM / SS) odgovara morfologiji zuba, sa nekoliko abrazivnih površina i kontaktnih površina. Baza zuba je smanjena.

Montažni sto 140
Sagitalni nagib, bočni kondilar (SKN) FH 30° - CE 15°
Bennetni ugao 10° ISS 1,0 mm
Laterotrusion LRT -20°
Retrusion RT 1,0 mm




Sekundarna morfologija abrazivna



Sekundarna morfološka abrazivna (SLA / SMA / SSA) odgovara morfologiji zuba koji predstavlja višestruke abrazivne površine i višestruke kontaktne površine. Zahvaljujući novom obliku tunela na bazi zuba, rad na korekciji se može smanjiti do 70 - 80%. U svim protetskim radovima može se primeniti ova abrazivna sekundarna morfologija.



Montažni sto 140
Sagittal nagib, bočni condilar (SKN) FH 25 ° - CE 10 °
Bennet kut 15 ° ISS 1,5 mm
Laterotrusion LRT -20 °
Retrusion RT 1,0 mm

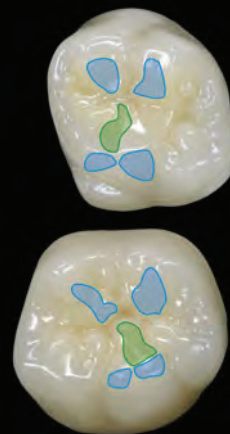

Tribos 501 c

Sekundarna morfologija abrazivna (kompozitna)



Tribos 501 c su:

- mnogo tvrdi od PMMA zuba (oko 30%)
- naročito otporni na abraziju
- otporni na pritisak i sa visokom otpornošću na lom

